



FDMNES

User's Guide

FDMNES

Yves Joly

yves.joly@neel.cnrs.fr
Institut Néel, CNRS, BP 166
38042 Grenoble Cedex 9, France

March 2019

Outline

Introduction	5
A) General Presentation	7
B) Main indata file	11
C) Convolution	53
D) Parameter optimization	63
E) Extraction of DAFS scan and spectra	67
F) Unit cell multiplication	69
G) FDMX user's guide	71
H) 2D diffraction	75
List of the <i>fdmnes</i> keywords	83

Introduction

The FDMNES program calculates the spectra of different spectroscopies related to the real or virtual absorption of x-ray in material. It gives the absorption cross sections around the ionization edges, that is in the XANES energy range. Calculations can be performed along all the conditions of linear or circular polarization. In the same way, it calculates the structure factors and intensities of anomalous or resonant diffraction spectra (DAFS or RXD) for 3D diffraction and for surface diffraction (SRXRD). FDMNES also allows the comparison of the simulated spectra to experimental ones with the help of objective criteria.

FDMNES is mainly a fully relativistic DFT-LSDA code. Optionally Hubbard correction (LSDA+U) can be used. It uses two techniques. The first one is based on the Finite Difference Method (FDM) to solve the Schrödinger equation. In that way the shape of the potential is free and in particular avoid the muffin-tin approximation. The second one uses the Green formalism (multiple scattering) on a muffin-tin potential. This approach can be less precise but is faster. The program includes also a multi-electronic extension using the “Time-Dependent DFT” with a local kernel. The program is symmetrized. Symmetry operations are calculated automatically.

The next section contains the information for practical purposes to run the program and in particular the description of the indata files. An introduction to x-ray spectroscopies is available in French.

The FDMNES program can be freely downloaded at the web address:

<http://www.neel.cnrs.fr/fdmnes>

In case of publication related to the use of the program thanks to cite:

O. Bunau and Y. Joly

Self-consistent aspects of x-ray absorption calculations

J. Phys. : Condens. Matter **21**, 345501 (2009).

The FDMNES program highly benefited from the scientific contribution of Calogero Natoli who has provided a constant and essential help. He is in particular at the origin of all the developments using the multiple scattering theory and the extensions to resonant diffraction and magnetism. Oana Bunau realized the extension towards TD-DFT and participated to the inclusion of self-consistency and Hubbard correction and now Yvonne Soldo-Olivier develops the surface x-ray diffraction. The program also benefited from the expertise of Delphine Cabaret, Hubert Renevier, Sergio Di Matteo, Christian Brouder, Aline Ramos, Stéphane Grenier, Emilio Lorenzo and Yvonne Gründer, without whom different advances would have not been realized. Finally, this work has been made greatly easier with the support of Denis Raoux.

From March 2015, Guda, Soldatov *et al.* (Rostov-on-Don, Russia) make the code far faster by the use of the MUMS library which solves systems of linear equations for sparse matrices. This method is now the default option used when running in the finite difference method mode. When using it thanks to cite also:

Sergey A. Guda, Alexander A. Guda, Mikhail A. Soldatov, Kirill A. Lomachenko ,
Aram L. Bugaev, Carlo Lamberti, Wojciech Gawelda, Christian Bressler, Grigory
Smolentsev, Alexander V. Soldatov, Yves Joly

*"Optimized Finite Difference Method for the Full-Potential XANES Simulations:
Application to Molecular Adsorption Geometries in MOFs and Metal-Ligand
Intersystem Crossing Transients"*

J. Chem. Theory Comput. **11**, 4512-4521 (2015).

and

P. R. Amestoy, A. Guermouche, J.-Y. L'Excellent and S. Pralet

Hybrid scheduling for the parallel solution of linear systems

Parallel Computing Vol 32 (2), pp 136-156 (2006).

From October 2015, the code also includes the FDMX extension from J. Bourke and
Ch. Chantler. When using it, thanks to cite also:

Jay Daniel Bourke, Christopher Thomas Chantler and Yves Joly

*"Extended X-ray Absorption Fine Structure Calculations Using the Finite Difference
Method"*

J. Synchrotron Rad. **23**, 551-559 (2016).

To have an overview of X-ray absorption spectroscopy, we recommend the reading of:
X-Ray Absorption and X-ray Emission Spectroscopy: Theory and Applications
Edited by J. A. van Bokhoven and C. Lamberti, Wiley (2016).
ISBN: 978-1-118-84423-6.

The work presented here concerns more specifically the chapter 4:
"Theory of X-ray Absorption Near Edge Structure"
by Yves Joly and Stéphane Grenier.

Another of our publications gives the main concepts of resonant diffraction:
"Basics of Resonant Elastic X-ray Scattering theory"
by S. Grenier and Y. Joly,
J. Phys.: Conference Series **519**, 012001 (2014).

The simulation of surface resonant x-ray diffraction and x-ray Raman spectroscopy are
respectively given in:

Yves Joly, Antoine Abisset, Aude Bailly, Maurizio De Santis, Farid Fettar, Stéphane Grenier,
Danny Mannix, Aline Y. Ramos, Marie-Claire Saint-Lager, Yvonne Soldo-Olivier, Jean-Marc
Tonnerre, Sergey A. Guda, and Yvonne Gründer.

"Simulation of surface resonant x-ray diffraction"

J. Chem. Theory Comput. **14**, 973-980 (2018).

DOI: 10.1021/acs.jctc.7b01032.

Yves Joly, Chiara Cavallari, Sergey A. Guda, and Christoph J. Sahle

"Full potential simulation of x-ray Raman scattering spectroscopy"

J. Chem. Theory Comput. **13**, 2172-2177 (2017).

DOI: 10.1021/acs.jctc.7b00203.

A- General presentation

I- Computer configuration

FDMNES run on all the computers, under LINUX or Windows. The programming language is Fortran 90. Executables are furnished for Windows and Linux 64 bits for sequential calculations.

Users can also compile the code themselves. The code needs then MUMPS, LAPACK and BLAS libraries. When they are not available, it is possible to use a gaussian solver routine furnished with the other routines

The code can be run on parallel processors if it is linked with the MPI library.

II- The package

Several groups of files can be downloaded: the *fdmnes* program itself and a set of examples of indata files.

- *fdmnes_win64.exe* is the executable for Windows 64 bits,
- *fdmnes_win32.exe* is the executable for Windows 32 bits,
- *fdmnes_linux64* is the executable for linux 64 bits,
- *fdmfile.txt* is an indata file,
- *prog* is a directory containing all the subroutines (*main.f90*, *general.f90*...), and the files included during compilation ("*mpif.h*"...).
- *Sim/Test_stand/in* is a directory containing a set of examples of *fdmnes* indata files.
- *Doc* is a directory contains information on the code, the manual and the the file « FDMNES_Modifications.txt » where are listed the modifications done in the code along the time.

An executable for the running in parallel can also be downloaded separately.

IV- Compilation

For those who do not want or cannot use the executables of the package, it is possible to compile the fortran routines. The code is more efficient when running with MUMPS library. All the furnished routines, but *mat_solve_gaussian.f90*, *sub_util.f* and *not_mpi.f90*, must then be compiled and linked with the call to this library and the ones whom it depends (SCOTCH, METIS, BLAS, LAPACK...). The files included in this step (*mpif.h*, ...) must be set in the same directory (or in a dedicated one).

When this library is not available, one exchanges *mat_solve_MUMPS.f90* by *mat_solve_gaussian.f90*, *sub_util.f* and *not_mpi.f90*, and one compiles and links without call to external library. Nevertheless by this way, calculations will be far more CPU consuming. *sub_util.f* can be replaced by BLAS and LAPACK.

Two examples of « makefile » for linux are provided.

V- Parallelization

Thanks to Sergey Guda, Keisuke Hatada, Kuniko Hayakawa and Rainer Wilcke, the users having the access to a cluster of computers can, using the MPI library, run the program in parallel mode. For this purpose, one has to delete the files "*mpif.h*" and "*not_mpi.f*" when compiling and makes the call to the corresponding library.

VI- Running

After compilation, the program can be run following the usual procedure available on your system.

As soon as the program is running, it calls the file "*fdmfile.txt*". This file must also be in the same directory than the executable file. It only contains the number of independent calculation to perform, followed the name of the indata file of each of these calculations. For example:

! Indata file for fdmnes

1

Sim/cu/in/cu_inp.txt

→ number of indata files

→ name of the indata file

When running in parallel, with MPI, one can use the executable given for this purpose, or another one compiled by the user with the MUMPS library. There are 2 levels of parallelization in FDMNES. The first one is a simple one inside the loop over energy. The second one, optional works inside the MUMPS library. It is useful only when working using the finite difference method.

Typically in the first case, the command to send the job is:

```
mpirun -np 16 fdmnes
```

The calculation uses then 16 processors in parallel.

To use the 2 parallelism one typically must write:

```
mpirun -np 16 -x HOST_NUM_FOR_MUMPS=2 fdmnes
```

In that case, 16 processors are used, 8 energies will be calculated in parallel. For each of them, 2 processors run the MUMPS library.

VII- Troubleshooting

When the program stops without reason, check if an "*fdmnes_error.txt*" file has been created. When created, this file contains a message explaining the trouble. Most often it is due to an error in the indata file. It can also be due to the fact that the indata file is not found. One has to check several points:


- 1) The called files are in the good directory. The name must contain all the paths.
- 2) The extensions are correct.
- 3) Their names are OK (under linux, upper and lower cases must be respected).

- 4) When using downloaded indata files, some problems of compatibility between systems can occur. It can be better to write again completely these files.

When the program stops without “*fdmnes_error.txt*” file, it can be due to a problem of space memory. Sometimes one gets a message with “stacking fault”. In this case try again putting in the indata file the keyword “Memory_save”. This keyword can be useful when there are many non-equivalent atoms. An approximation (in fact very good) is then done on the potential calculation. This option saves some memory space.

VIII- Structure of the calculation

For a complete calculation one has the following scheme:

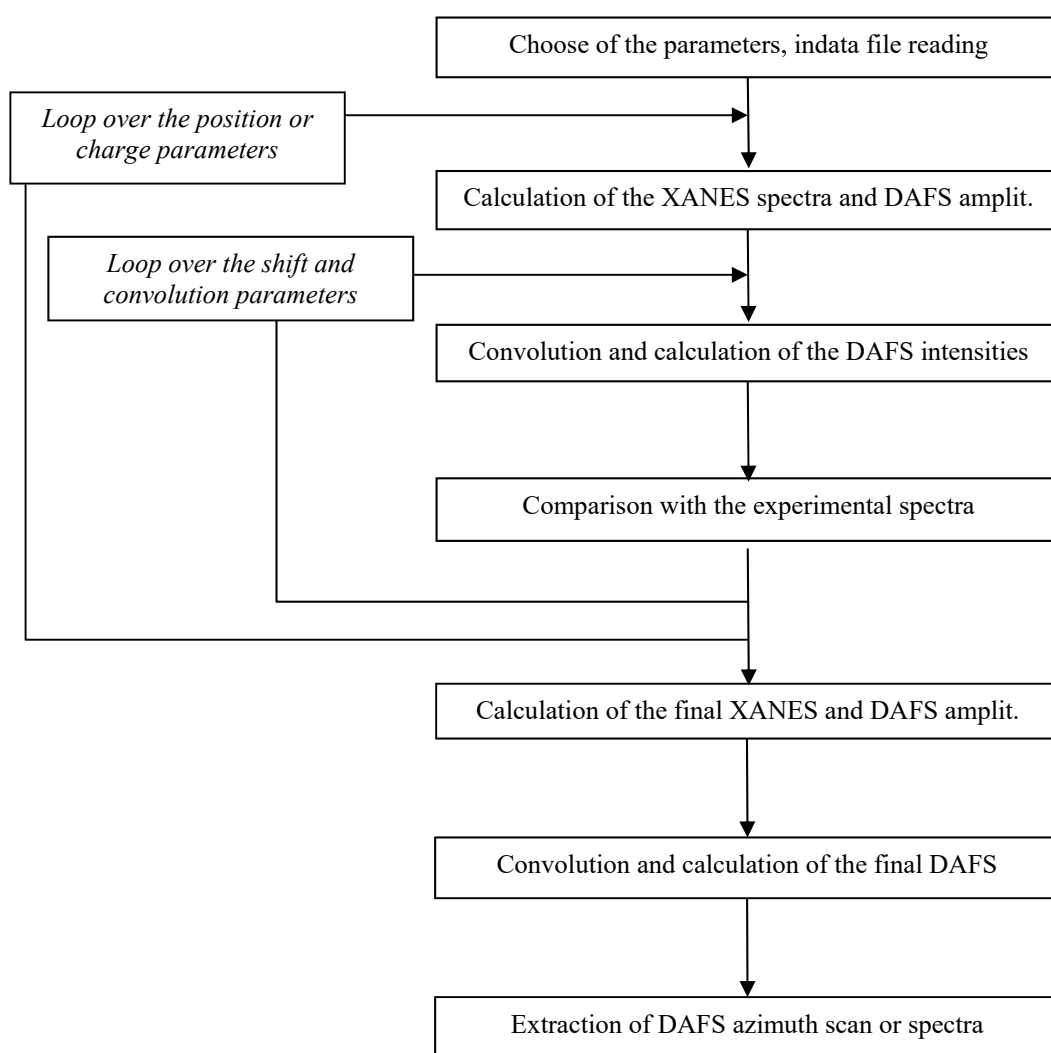

The program allows the calculation of spectra from a grid of parameters (position, charge...), to convolute them with other parameters (convolution width, energy shift...), then to compare them to experimental spectra with the help of objective criteria. These different steps can be performed together or separately. The comparison with the experimental spectra is also not mandatory. When used, the parameter fit must be performed with care. In practice

many calculations are limited to the step “XANES and DAFS calculation” and “Convolution and calculation of DAFS intensities”. These two steps can also be performed together or separately.

In the output files, the absorption cross section are in Mbarn ($1 \text{ Mbarn} = 10^{-18} \text{ cm}^2$) and summed up over the atoms of same chemical specie in the unit cell or in the cluster. To convert in number of electron one has to multiply by $C = \frac{\hbar\omega_{eV}}{800\pi^2 a_0^2 \alpha R} = 0.004555352 \times \hbar\omega_{eV}$, where R , α and a_0 are respectively the Rydberg constant, the fine structure constant and the Bohr radius in Angstrom, $\hbar\omega_{eV}$ is the photon energy in eV. One has also to divide by the number of atoms if one wants the result per atom. The intensities of the reflections are in square of number of electron.

The next chapter treats about the principal indata file for the step “XANES and DAFS calculation”. Generally this file is sufficient to describe all the necessary data for the calculation because the program calculates its atomic bases and the potential. Nevertheless, the user can prefer use its own atomic bases or uses directly the potential calculated by the band structure program FLAPW WIEN-2k. In both cases, some other files must be furnished. They are described further in the manual. The indata necessary for the steps “Convolution”, “comparison with the experimental spectra” and “Extraction of azimuth scan or spectra” can be set in the same indata file, but they are explained separately in the chapter C, D and E.

B- Main indata file

I- General Structure

It contains most of the inputs necessary for the calculation. All the data in input and output are in Angstrom and electron-Volt. Many options are chosen by default. One can modify or add other options using keywords. Text can be in upper or lower case. Blank lines or beginning by “!” are not taken into account. Between number, one must put at least one blank. When getting problem when opening these indata files, one has to check if their name is correct. Moreover some compilers do not like files written under other system (MAC, DOS, LINUX...). In case of difficulties when the program wants to open one of these downloaded files, it can be useful to completely write them again.

The indata file contains several blocks of data, each ones starting with a specific keyword. The end of the indata file is noted by the keyword "*End*". Whatever is after is not read. Here comes an example of indata file:

! Fdmnes indata file

! Calculation for the copper K-edge in copper cfc

Filout

Sim/Cu/Cu_out

→ Name of the output files (without extension)

Range

-2. 0.2 5. 0.5 10. 1. 40.

Radius

3.0

→ Cluster radius

Crystal

3.610 3.610 3.610 90. 90. 90.

29 0.0 0.0 0.0

29 0.5 0.5 0.0

29 0.5 0.0 0.5

29 0.0 0.5 0.5

→ crystal structure

→ Unit cell parameters (Å and degrees) : a, b, c, α , β , γ

→ Atomic number, position

Convolution

→ to get convoluted spectra

End

→ end of the indata file

Two blocks are necessary for any calculations. The first one starts with the keyword "*Radius*" followed by the value of the radius of the sphere surrounding the area of calculation. The second one is necessary to describe the material structure. When that one comes from the output of the WIEN-2k package, it begins with the keyword "*FLAPW*". In the other cases, the molecule structure or the elementary unit cell, in case of a 3D structure, is described with the keyword "*Molecule*" or "*Crystal*", followed by the unit cell parameter and atom positions.

All the keywords related to the convolution or to the fit of the parameters are treated in chapter C and D.

Output file names

By default the output file name is *fdmnes_out*. This name can be modified by the use of the keyword "*filout*" followed by the name we want (without extension). Then one gets several output files with the extensions:

<i>_bav.txt</i>	output file giving details
<i>.txt</i>	contains only the spectra by column
<i>_nrixs.txt</i>	contains only the spectra by column for NRIXS simulations

If a calculation is performed on several non-equivalent crystallographic sites, one gets the extensions:

_i.txt, _j.txt ... in which *i* and *j* are the index of the sites (see keyword *absorber*)

In option or depending on the type of calculation, one can also get the files:

<i>_conv.txt</i>	convoluted spectra scan (keyword <i>Convolution</i>).
<i>_scan.txt</i>	dafs versus angles for azimuthal scan (keyword <i>DAFS</i>).
<i>_sda.txt</i>	state density for the atom number <i>a</i> (keyword <i>Density</i>).
<i>_atoma.txt</i>	results for one atom at position number ' <i>a</i> ' (keyword <i>Allsite</i>).
<i>_atoma_scan.txt</i>	<i>DAFS</i> scan results for the atom <i>a</i> (keyword <i>Allsite</i> and <i>DAFS</i>).
<i>_tddft.txt</i>	output with the TDDFT option (keyword <i>Tddft</i>).
<i>_tddft_scan.txt</i>	azimuthal scan in the TDDFT option (keyword <i>DAFS</i> and <i>Tddft</i>).
<i>_tddft_conv.txt</i>	convoluted spectra in TDDFT (keyword <i>Convolution</i> and <i>Tddft</i>).

For the analysis of the spherical tensors (keyword *Spherical* and *Sphere_all*)

<i>_sph_atoma.txt</i>	spherical tensors of the atom <i>a</i> .
<i>_sph_atoma_int.txt</i>	integral of the spherical tensors of the atom <i>a</i> .
<i>_sph_signal_atoma_xan.txt</i>	contribution of each atomic spherical tensor on the average xanes signal.
<i>_sph_signal_atoma_poli.txt</i>	contribution of each atomic spherical tensor on the xanes polarisation number <i>i</i> .
<i>_sph_signal_atoma_rxsi.txt</i>	contribution of each atomic spherical tensor on the <i>DAFS</i> reflection number <i>i</i> .
<i>_sph_signal_xtal_xan.txt</i>	contribution of the crystal spherical tensor on the average xanes signal.
<i>_sph_signal_xtal_rxsi.txt</i>	contribution of the crystal spherical tensor on the <i>DAFS</i> reflexion number <i>i</i> .
<i>_sph_xtal.txt</i>	spherical tensors of the crystal.
<i>_sph_xtal_int.txt</i>	integral of the spherical tensors of the crystal.
<i>_sph_xtal_rxsi.txt</i>	spherical tensor of the crystal for the <i>DAFS</i> reflexion number <i>i</i> .

For the analysis of the cartesian tensors (keyword *cartesian*):

<i>_car_atoma.txt</i>	cartesian tensors of the atom <i>a</i> .
<i>_car_xtal.txt</i>	cartesian tensors for the crystal
<i>_car_xtal_rxsi.txt</i>	cartesian tensors for the crystal for the <i>DAFS</i> reflection number <i>i</i> .

II- Basic keywords

1) Output file names

The different output files have names with the same root. The extensions automatically added depending on the chosen option. To define this root use :

Filout

Sim/Cu/Cu_out

→ Name of the output files (without extension)

The files can eventually be in a subdirectory.

2) Radius of the cluster

The final states are calculated inside a sphere, whose radius is defined with the keyword "*Radius*". Only the atoms inside this sphere are considered. By default, the sphere is centered on the absorbing atom.

Radius

→ Mandatory keyword preceding the radius of the cluster.

3.5

→ value in Angstrom of the cluster radius.

For calculations on large energy range (typically used in EXAFS), it can be convenient using a decreasing radius when energy increases. For this purpose, one can set energy ranges with decreasing radius:

Radius

8.3 100 7 250 5 500 4 → 8.3: radius up to 100 eV, 7 radius between 100 and 250 eV...
4 : radius beyond 500 eV

3) Cluster or crystal structure

Under "*Crystal*" or "*Molecule*" stand all the data describing respectively the unit cell or the molecule. If the calculation is done using the FLAPW output, this block is useless because the structure is read in one of the FLAPW output files. Under the keyword, come the mesh parameters (Å) and the angles (degrees). Then come all the atoms (and not only the non-equivalent ones, but when is specified the space group using the keyword "*Spgroup*"). By default and in the absence of the keyword "*Absorber*", the absorbing atom chemical specie corresponds to the first atom in the list.

Example 1: cfc copper crystal:

Crystal

```
3.610 3.610 3.610 90. 90. 90.
29 0.0 0.0 0.0
29 0.5 0.5 0.0
29 0.5 0.0 0.5
29 0.0 0.5 0.5
```

→ Crystal structure

→ a, b, c, α, β, γ

→ Atomic number, position

Example 2: FeO₆ octahedron:

Molecule

	1.900	1.900	1.900	90.	90.	90.	→ a, b, c, α , β , γ
26	0.0	0.0	0.0				→ Atomic number, position
8	1.0	0.0	0.0				
8	-1.0	0.0	0.0				
8	0.0	1.0	0.0				
8	0.0	-1.0	0.0				
8	0.0	0.0	1.0				
8	0.0	0.0	-1.0				

The atomic structure can also be given in cylindrical or spherical coordinates. To use cylindrical coordinates, it is sufficient to give only two numbers under "*Molecule*". The program will understand they are *a* and *c* and that the positions of the atoms are given by *r*, ϕ and *z*. Thus to describe the same octahedron as previously:

Molecule

	1.900	1.900		→ a, c
26	0.0	0.0	0.0	→ Atomic number, position (<i>r</i> , ϕ and <i>z</i>)
8	1.0	0.0	0.0	
8	1.0	180.0	0.0	
8	1.0	90.0	0.0	
8	1.0	-90.0	0.0	
8	0.0	0.0	1.0	
8	0.0	0.0	-1.0	

To use spherical coordinates, only one number (*a*) must be set after "*Molecule*". The position of the atoms is then given by *r*, θ , ϕ . For the same octahedron:

Molecule

	1.900			→ a
26	0.0	0.0	0.0	→ Atomic number, position (<i>r</i> , θ , ϕ)
8	1.0	90.0	0.0	
8	1.0	90.0	90.0	
8	1.0	90.0	180.0	
8	1.0	90.0	270.0	
8	1.0	0.0	0.0	
8	1.0	180.0	0.0	

When working with 2D film, one just has to replace « *Crystal* » by « *Film* ».

It is possible, in the non-magnetic case, to specify only the non-equivalent atoms. Then one has to give the space group under the keyword "*Spgroup*". The complete name, as in the international table, must be given. For example, for magnetite, one gets:

Spgroup

Fd-3m:1 → it is also possible to write 227:1

Crystal

```
8.3940 8.3940 8.3940 90.0 90.0 90.0
```

```
26 .6250 .6250 .6250 ! Fe 16d
26 .0000 .0000 .0000 ! Fe 8a
8 .3800 .3800 .3800 ! O 32e
```

Note that for the calculation of the symmetry, it is important to define the atom position with a sufficient number of digits (say 10). For instance, for graphite one has to write:

Spgroup

```
P63mc
```

Crystal

```
2.456 2.456 6.696 90. 90. 120. = a, b, c, alpha, beta, gamma
6 0.0 0.0 0.0
6 0.3333333333 0.6666666667 0.0
```

Writing for example 0.3333 in the last line, would create false atoms. Program would stop immediately with an error message.

It is also possible to have an occupancy rate not full for the atoms. For this one must write the keyword "Occupancy" and a fifth column after the atom position containing this weight:

Occupancy

Crystal → or **Molecule** or **Film** or **Surface** or **Interface**

```
2.456 2.456 6.696 90. 90. 120.
6 0.0 0.0 0.0 0.8 → weight 80 %
6 0.25 0.25 0.0 1. → weight 100%
```

During the building of the cluster, when 2 atoms are at the same position, only the one with the highest weight is kept. The final signal is taken with the given weights.

It is possible to use directly the structure given in *pdb* or *cif* files. These ones contain all the geometry, the space group... In this case use, one of the following way:

Pdb_file (or **Film_Pdb_file**, when working with a 2D film)

```
File_name.pdb
```

Cif_file (or **Film_Cif_file**, when working with a 2D film)

```
File_name.cif
```

It is also possible to describe a surface structure with a 2D periodicity using the keywords « *Surface* », « *Film* » or « *Bulk* ». See these keywords presented in section "2D Diffraction". They can also be used for XANES.

4) Atomic electronic densities

An electronic configuration is used by default for all the atoms. It is possible to modify it by the use of the keyword "*atom*". Moreover under "*crystal*" or "*molecule*" one must not anymore put the atomic number, by the atom type number. For example in case of a FeO₆ octahedron with the Fe 3d⁶4s² and O 2s²2p⁴ configuration:

Atom → keyword preceding the atomic electronic densities
 26 2 3 2 6. 4 0 2. → atomic number of the chemical specie of type 1, number
 8 2 2 0 2. 2 1 4. of valence orbital and (n,l,po) of each of these orbitals

molecule

	1.900	1.900	1.900	90.	90.	90.	→ a, b, c, α, β, γ
1	0.0	0.0	0.0				→ Atom type, position
2	1.0	0.0	0.0				
2	-1.0	0.0	0.0				
2	0.0	1.0	0.0				
2	0.0	-1.0	0.0				
2	0.0	0.0	1.0				
2	0.0	0.0	-1.0				

Important remark: contrary to what one can think, the formal charges attributed to the atoms in the ionic compounds are far from the true charge. Thus one has to perform exchange of charge between atoms with care and in a moderate way. A good technique is, for example for 3 d elements, the good number of "d" electron, following the formal charge, but keeping the neutral atom, putting electrons in the large radius 4s or 4p orbitals.

Another possibility permits to keep the atomic number in front of each atom in the list after *Crystal* or *Molecule* and to give the electronic configuration only for some of them. For this use the keyword:

Atom_conf

3 1 2 2 3 2 5. 4 1 1. → nbr of atom with the configuration then index of these atoms,
 1 3 2 3 2 6. 4 1 0. Then number of orbital and configuration of them

Crystal

4. 4. 4. 90 90 90
 26 0. 0. 0.
 26 0.5 0. 0.
 26 0. 0.5 0.
 26 0. 0. 0.5

In the example above the atoms 1 and 2 of the list under « *Crystal* » have the configuration 3d⁵4p¹. The atom 3 has the configuration 3d⁶. The remaining atoms have the default configuration.

When one wants to give a configuration for a doping element (see keyword « *Doping* »), one must write « 0 » for the atom index:

Atom_conf

1 0 2 3 2 5. 4 1 1. → nbr of atom (1), then index = 0

5) Absorbing atoms

All the atoms present in the structure participate to the absorption or scattering. By default, the calculated spectra correspond to the sum of the scattering produced by all the atoms of the same atomic number than the first one in the list under "*Crystal*" or "*Molecule*".

When one wants to select some of the absorbing atoms, one has to use the keyword "*Absorber*". For example, if one wants that the absorbing atom is the n^{th} of the list and only this one (under "*Crystal*" or "*Molecule*"), put:

Absorber

3 → absorbing atom number (here the 3rd in the list).

If this example goes with the FeO₆ structure defined just above, the result in the output file would correspond to a calculation of an oxygen atom summed over the 6 atoms because they are equivalent by symmetry. If one wants the result from a unique site, use the keyword "*Symsite*" or "*Allsite*" (see further on).

If there are non-equivalent sites, one has to perform independent calculations. For this, one must specify each atom number under "*Absorber*". The name of each of the corresponding output file get the suffix "*_n*" in which n is the atom number.

Filout

copper_out

Absorber

1 5 → atom numbers whom results will in output file *copper_out_1* and *copper_out_5*

The summation over these different sites and the eventual energy shifts is to be performed during a following step.

Another possibility is to define the atomic numbers of the absorbing atoms. We can have 1 or several atomic numbers (but of the same edge, K...) This is done with the keyword:

Z_absorber

26 → all the atoms with $Z = 26$ are absorbing atoms

Z_absorber

26 27 → all the atoms with $Z = 26$ or $Z = 27$ are absorbing atoms

6) Energy range

The energy range E that one defines in the *indata* is the energy of the photoelectron relative to the Fermi level.

By default the energy range is -5 to 60 eV by 0.5 eV step. One can change the range, the step or even have a variable step using:

Range → keyword for the energy range
1. 0.5 60. → E_{\min} , step, E_{\max}

Other example with variable step:

Range
1. 0.1 10. 0.5 20. 1. 60.00 → E_{\min} , step, $E_{\text{intermediate}}$, step ...

To get a continuously increasing step (k step constant) put:

Rangel
1. 0.1 200. → E_{\min} , step at the Fermi level, E_{\max}

By default, the output energy range is relatively to the Fermi level. If one wants that the output energy is the photon energy put the keyword:

Energpho

7) Multiple scattering mode

If one wants to calculate in the multiple scattering mode use the keyword:

Green

Then the potential is automatically a muffin-tin one. The mode is faster than the finite difference method, so one has to use it first.

8) Threshold type

By default the threshold is the K one. For other threshold put the keyword:

Edge → keyword preceding the threshold type
L1 → threshold (K, L1, L2, L3, M1 ...)

It is possible in a single run to calculate 2 edges with the same initial (n,l), that is the edges L2 and L3 or M2 and M3 or M4 and M5 ... For this write:

Edge → keyword preceding the threshold type
L23 → threshold L2 and L3, or M23, M45, N23, N45.

9) Multipolar expansion

By default only the transition electric dipolar component (E1E1) is calculated. This is modified by the keywords:

Quadrupole → quadrupolar calculation (E1E2 and E2E2)

<i>Octupole</i>	→ octupolar calculation (E1E3 and E3E3)
<i>Dipmag</i>	→ magnetic dipole calculation (E1M1) and (M1M1)
<i>E1E2</i>	→ calculation of E1E2
<i>E1E3</i>	→ calculation of E1E3
<i>E2E2</i>	→ calculation of E2E2
<i>E3E3</i>	→ calculation of E3E3
<i>E1M1</i>	→ calculation of E1M1
<i>M1M1</i>	→ calculation of M1M1
<i>No_E1E1</i>	→ No calculation of the dipolar-dipolar component (E1E1)
<i>No_E2E2</i>	→ No calculation of the quadrupolar-quadrupolar component
<i>No_E1E2</i>	→ No calculation of the dipolar-quadrupolar component
<i>No_E1E3</i>	→ No calculation of the dipolar-octupolar component
<i>Dip_rel</i>	→ To take into account the spin-orbit transition term in E1E1

10) Polarization and dichroism

By default the calculation is performed along 1, 2 or 3 orthogonal polarizations in the dipole mode (and up to six in quadrupolar) depending of the symmetry, and only the average corresponding to a powder is given. To get the spectra for the different polarization, use:

Polarize

To get linear or circular dichroism, one must specifies polarization and wave vector (in quadrupolar) orientations. These vectors must then be given in the non-orthogonal, but normalized crystal basis:

Polarize

```
1.0 1.0 0.0
1.0 -1.0 0.0
0.0 0.0 1.0
```

For a quadrupolar calculation, one has to specify the wave vector after the polarization:

Polarize

```
1.0 1.0 0.0 0.0 0.0 1.0
1.0 -1.0 0.0 0.0 0.0 1.0
0.0 0.0 1.0 0.0 0.0 0.0
```

Each line contains the polarization vector, the wave vector. If the wave vector is zero this polarization is calculated in the dipolar approximation.

It is possible to perform an average on different polarizations adding a new number at the end of the line is the corresponding weight. If at least two weights are non zero, in the output, there will be a new column with the weighted average of the different polarization.

Polarize

```
1.0 1.0 0.0 0.0 0.0 1.0 1.
1.0 -1.0 0.0 0.0 0.0 1.0 1.
0.0 0.0 1.0 0.0 0.0 0.0 0.
```

If one wants to have a circular polarization, one just has to put this one to zero followed by the value of the wave vector:

Polarize

0.0 0.0 0.0 0.0 0.0 1.0

In the output one will have the calculation in right polarization, then the left polarization, then the difference. In this example, the calculation corresponds to the polarization $x+iy$ and $x-iy$.

To calculate the transmission through a sample of a specific thickness and to check change in polarization to the sample, it is necessary to calculate an absorption σ - π matrix. This is done with the keyword:

Mat_polar

0. 0. 1. 1. 0. 0. → polarization σ and wave vector

See also the keywords: “*Sample_thickness*”, “*Stokes*” and “*Stokes_name*” to define the sample thickness, incoming polarization matrix and analyzer angles.

II-10) Doping

It is possible, more specifically for crystal to calculate the absorption of a doping element set in substitution in place of another one. For this one uses the keyword « Doping », followed by the atomic number of the doping element and the index of the atom to substitute in the list under « *Crystal* »:

Doping

27 1 → substitution by Cobalt of the atom in the first site.

The doping atom is supposed to be at low concentration, thus the clust built around it is the same than the one given by the crystal. Symmetries are kept.

II-11) Getting back the electronic structure from a previous calculation

It is possible to get back the electronic structure of the absorbing atoms from a previous calculation (it is in the *fdmnes_out_bav.txt* type file) and to compute other polarizations, or ask for the density of state (keyword Density), or make a TDDFT calculation different or not done in the previous calculation, putting the key word “*Extract*” followed by the corresponding name:

Extract

Sim/calcul_prec_bav.txt

III- Non resonant X-ray inelastic scattering (NRIXS)

The non-resonant inelastic scattering (NRIXS or X-Raman) is very close to XANES. The first order term is equivalent to the dipolar approximation in XANES. One can obtain a more precise simulation, depending on the momentum transfer q , using the keyword:

NRIXS

3. 5.5 9. → the different values of the modulus of q in \AA^{-1}

The corresponding output file will have the extension "_nrixs".

When one wants to simulate an experiment performed on a mono-crystal one must specify also the direction of vector q . This is done replacing the keyword "*NRIXS*" by "*NRIXS_mono*", followed line by line by the modulus of q and the direction of q which is expressed in the unit cell basis. This direction does not need to be normalized:

NRIXS_mono

3. 1. 0. 0. → modulus of q , then direction of q

3. 1. 1. 1.

3. 1. 2. 2.5

3. → when the direction vector is not given (or equal to (0,0,0)) this

5. 1. 0. 0. corresponds to the calculation of powder, like when using « NRIXS »

5. 1. 1. 1.

If one wants the results resolved in ℓ , the Bessel quantum number associated to the expansion of $\exp(iq \cdot r)$ in the transition operator, write:

All_nrixs

By default the maximum value of ℓ is 2. To change it, write :

Lmax_nrixs

3 → ℓ maximum value.

IV- Anomalous or resonant diffraction

1) Reflections and polarization

In case of anomalous (or resonant) diffraction (DAFS, DANES, RXS, RXD) calculation, put the keyword "*DAFS*" (or "*RXS*") followed by the index of the beams to calculate. The orientation of the polarization and wave vector can be described by different way. When working in σ - σ , σ - π , π - π or π -s, or in circular polarization, put the number 1, 2, 3, 4 or 5 for the polarization respectively σ , π , circular right, circular left or linear along a general direction, in input then in output.

DAFS

0 0 2 1 2 45. \rightarrow reflection indices, σ , π , azimuth
 0 0 2 1 1 45. \rightarrow reflection indices, σ , σ , azimuth

When the polarization is linear but not σ or π , but with an angle α such that $\alpha = 0^\circ$ when it is σ and 90° when it is π , one must write:

DAFS

0 0 2 1 0. 5 10. 45. \rightarrow reflection indices, σ , angle, rectilinear, angle, azimuth
 0 0 2 1 0. 5 -10 45.

Note that in this case, one must specify both incoming and outgoing polarization angles, even when one of them is σ , π or circular. When it is circular, the angle is not taken into account.

Let's take \mathbf{a} , \mathbf{b} and \mathbf{c} , the unit cell vectors, $\mathbf{i}_a = \mathbf{a}/|\mathbf{a}|$, $\mathbf{i}_c = \mathbf{c}/|\mathbf{c}|$, and \mathbf{Q} the normalized diffraction vector. We define the azimuthal angle, φ , using the base $(\mathbf{I}, \mathbf{J}, \mathbf{Q})$ such that:

$$\mathbf{I} = \mathbf{Q} \times \mathbf{i}_c \times \mathbf{Q}/|\mathbf{i}_c \times \mathbf{Q}|, \quad \mathbf{J} = \mathbf{Q} \times \mathbf{I}$$

When \mathbf{Q} is along \mathbf{c} , \mathbf{i}_c is substituted by \mathbf{i}_a . One then gets, versus the Bragg angle θ_B the incoming and outgoing wave vectors:


$$\begin{aligned} \vec{k}_i &= \cos \theta_B \cos \varphi \vec{I} - \cos \theta_B \sin \varphi \vec{J} - \sin \theta_B \vec{Q}, \\ \vec{k}_s &= \cos \theta_B \cos \varphi \vec{I} - \cos \theta_B \sin \varphi \vec{J} + \sin \theta_B \vec{Q}. \end{aligned}$$

For the σ and π polarizations, one gets:

$$\begin{aligned} \vec{\varepsilon}_\sigma &= \sin \varphi \vec{I} + \cos \varphi \vec{J} \\ \vec{\varepsilon}_\pi &= \vec{k}_i \times \vec{\varepsilon}_\sigma = \sin \theta_B \cos \varphi \vec{I} - \sin \theta_B \sin \varphi \vec{J} + \cos \theta_B \vec{Q} \\ \vec{\varepsilon}_\pi &= \vec{k}_s \times \vec{\varepsilon}_\sigma = -\sin \theta_B \cos \varphi \vec{I} + \sin \theta_B \sin \varphi \vec{J} + \cos \theta_B \vec{Q} \end{aligned}$$

By this way, φ increases when the sample is rotated trigonometrically. For example with an orthogonal unit cell: reflection (h,0,0), $(\mathbf{I}, \mathbf{J}, \mathbf{Q}) = (\mathbf{i}_c, -\mathbf{i}_b, \mathbf{i}_a)$

1 1 0. \rightarrow corresponds to σ - σ with polarization along $-\mathbf{i}_b$.
 1 1 90. \rightarrow corresponds to σ - σ with polarization along \mathbf{i}_c .
 reflexion (0,0, ℓ), $(\mathbf{I}, \mathbf{J}, \mathbf{Q}) = (\mathbf{i}_a, \mathbf{i}_b, \mathbf{i}_c)$
 1 1 0. \rightarrow corresponds to σ - σ with polarization along \mathbf{i}_b .
 1 1 90. \rightarrow corresponds to σ - σ with polarization along \mathbf{i}_a .

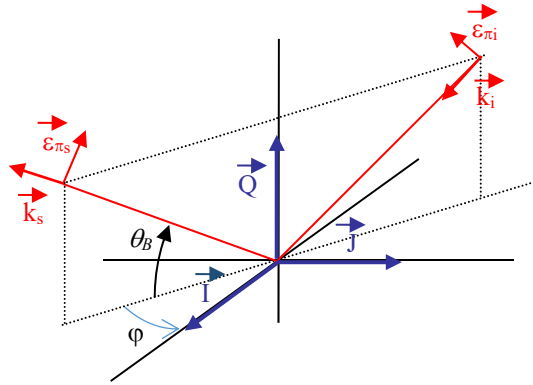

If one does not put the angle, this implies that one performs a φ scan and all the amplitudes are calculated for all the angles by 2° step. There is then a supplementary output files with the extension "_scan.txt".

It is also possible to perform a 360° scan by 2° step for the tilinear incoming or outgoing polarizations. For this, one has to write 10 in place of the σ - π notification:

DAFS

0 0 2 10 1 45. \rightarrow incoming polarization is scanned
 0 0 2 2 10 45. \rightarrow outgoing polarization is scanned

The first value (angle = 0°) corresponds to the σ polarization. 90° corresponds to π polarization.

It is possible to write the exact polarization directions. This can be useful for peculiar uses as in photoemission. In case of linear polarization write:

DAFS

0 0 0 \rightarrow Reflection index
 0. 0. 1. 0. 1. 0 $\rightarrow \epsilon_i, k_i$
 0. 0. 1. 0.7071 0.7071 0. $\rightarrow \epsilon_s, k_s$

For circular polarization, it is complex:

DAFS

0 0 0 \rightarrow reflection index
 0.7071 0.7071 0. 0. 0. 0. 0. 1. 0. $\rightarrow \epsilon_i, k_i$ ($\epsilon_{xr} \epsilon_{xi} \epsilon_{yr} \epsilon_{yi} \epsilon_{zr} \epsilon_{zi} k_{ix} k_{iy} k_{iz}$)
 0.7071 -0.7071 0. 0. 0. 0. 0. 1. 0. $\rightarrow \epsilon_s, k_s$

When one does not want phase term between atoms, for example to simulate photoemission, just write 0 0 0 as reflection index.

Note that it is possible to choose any other origin for the azimuth, just using another vector than \mathbf{i}_c to define the basis vector \mathbf{I} and \mathbf{J} . For this just write:

Zero_azim

0. 1. 1. \rightarrow vector in the direct crystal base

It is now possible to calculate the self-absorption corresponding to the incoming and outgoing photon polarizations. For this put the keyword:

In some cases, as in biology, one has a very large number of diffracted beams. It is then simpler to specify the beams to calculate using experimental files. In this case one must give the corresponding file names containing the experimental spectra, the orientation matrix giving the incoming polarization versus the unit cell and for each file the rotation angle. One thus has to write:

DAFS_exp

11.231547	3.519121	-71.823997	→ matrix
12.573563	-170.027649	-6.364535	
-86.558281	-5.978752	-11.279972	
0.			→ angle for the first file
MoFe/Refl_R5to2_1.txt			→ name of the first file
10.			→ angle for the second file
MoFe/Refl_R5to2_2.txt			→ name of the second file ...

Each experimental file contains some number of spectra in column along the format :

-1 16 28 -1 20 18	→ (h,k,l) for each reflections (integer)
Energy D[-1,16,28] D[-1,20,18]	→ one line of comment
7080 -3.19 -7.96	→ Energy, intensity for each reflections
7100 9.75 16.65	

One must have at least one space between each number, but no tabulation or coma between the numbers. The number of digit or space between the numbers is free.

For a calculation with *DAFS_exp*, by default one calculates and compares to the data, the value : $\sqrt{I_{\sigma\sigma}(Q) + I_{\sigma\pi}(Q)} - \sqrt{I_{\sigma\sigma}(-Q) + I_{\sigma\pi}(-Q)}$. If one wants to calculate and compare other values related to the intensity, write «*DAFS_exp_type*» with below a number between 1 and 4:

DAFS_exp_type

1	→ $I_{\sigma\sigma}(Q) + I_{\sigma\pi}(Q)$
2	→ $\sqrt{I_{\sigma\sigma}(Q) + I_{\sigma\pi}(Q)}$
3	→ $I_{\sigma\sigma}(Q) + I_{\sigma\pi}(Q) - I_{\sigma\sigma}(-Q) - I_{\sigma\pi}(-Q)$
4	→ default case : $\sqrt{I_{\sigma\sigma}(Q) + I_{\sigma\pi}(Q)} - \sqrt{I_{\sigma\sigma}(-Q) + I_{\sigma\pi}(-Q)}$

Self_abs

Then one gets in the output files, after each reflection two new columns containing the linear absorption coefficients, the unit being the μm^{-1} . They contain after the convolution also the absorption coming from the other atoms and the other edges. That is that the absorption before the edge is not zero. The new data allows the correction due to the self-absorption in order to compare with experimental spectra. The corrected spectra are, after convolution in columns following the non-corrected spectra. Their names begin with "Ic".

One can make a more sophisticated correction taking into account the birefringence. For this put the keyword:

Full_self_abs

It is then necessary to calculate all the reflections along the 4 polarizations $\sigma\sigma$, $\sigma\pi$, $\pi\sigma$ and $\pi\pi$ and in this order:

DAFS

2 0 0 1 1 0.

2 0 0 1 2 0.

2 0 0 2 1 0.

2 0 0 2 2 0.

In first approximation, this correction is equivalent to *Self_abs*. It allows the taking into account of the rotation of the polarization when the electromagnetic wave propagates into the media.

Note that from this, it is possible to calculate circular polarizations reflection intensity for any value of the Stokes parameters. This part is performed in the “convolution” step. See thus the corresponding paragraph.

By default the scan step is 2°. To modify it, use:

Step_azim

0.5 → step indegres of the azimuthal scan

2) Non resonant magnetic scattering

The non-resonant magnetic scattering is taken into account for the *DAFS*. This term can be decreased, or put to zero, using a multiplicative factor:

No_res_mag

0.8 → factor

An additional factor can be applied to the contribution from the orbital moment. This factor corresponds to $L/2S$ (and not L/S):

No_res_mom

-0.3 → factor for the orbital moment

By default, this factor is calculated using the Hund rules and multiplied by 0.2.

3) Temperature

It is possible to take into account the thermic disorder using the Debye model. For this purpose use the keyword :

Atom_B_iso

After keyword « Crystal », one then must add a fifth column after the atom positions the Debye parameter ($8\pi^2\langle u^2 \rangle$) in Å². If no number are given, the default zero value is given. When using also « *Occupancy* » keyword, the fifth and sixth columns contain the 2 corresponding parameters in the order given by the order of the keywords in the indata file.

For the same purpose but using as Debye parameter $\langle u^2 \rangle$ instead of $8\pi^2\langle u^2 \rangle$ use:

Atom_U_iso

4) Occupancy

In the same way, one can use an occupancy rate on each sites. For this purpose write:

Occupancy

Then after keyword “Crystal” in the fifth column, after the atom position, the occupancy rates, between 0 and 1 are given. When no number value, the default 1 value is used. When using also “*Atom_B_iso*” or “*Atom_U_iso*”, see just above for the order.

5) Common keywords with SRXRD

« *DAFS_2D* » allows the definition of reflections and polarizations using the operation modes corresponding to the diffractometer geometry with up to 6 circles. It is described in chapter SRXRD.

« *No_analyzer* » allows the calculation of reflections without analyzer. See also the chapter SRXRD.

V- Potential, magnetism and calculation technique

1) Relativistic calculation

By default, simulations are non-relativistic, without spin-orbit and non-magnetic for the valence and conduction band states (but relativistic with spin-orbit for core states). For heavy atoms (say $Z > 36$), in the structure, including when they are not absorbing atoms, it is recommended to make a relativist calculation using the keyword:

Relativism

To be more precise, it can be useful, to start from the Dirac equation, which contains the spin-orbit interaction. Such a calculation, fully relativistic can be done with the keyword:

Spinorbit

Note that then the simulations are typically 4 to 8 times longer and need 2 times more memory space. Beside the relativistic aspect, when scattering or absorption depend on conduction band electron spin-orbit, this keyword is mandatory. A spin-polarized simulation can then be done using this keyword plus the keyword "*Atom_conf*" or "*Atom*" (see further on).

With the keyword "Spinorbit", simulation is automatically relativistic. For a non-relativistic simulation without spin-orbit use the keyword:

Nonrelat

2) Spin polarized calculation

To make a spin polarized calculation, one can use the "*Spinorbit*" keyword, see above, or more simply, when spin-orbit can be neglected (what can be the case of L₂₃ edges but never of the K edges) put the keyword:

Magnetism

If the polarized potential comes from LAPW, "*spinorbit*" or "*magnetism*" must be before "*flapw*". When the calculation does not use the Wien-2k output, the use of the keywords "*Atom*" or "*Atom_conf*" is mandatory to specify the electronic configuration, different for the spin up and spin down parts. The orbital occupancy must be given for each spin in doubling the corresponding columns.

Example of metal fcc nickel:

Magnetism

Atom

```
28 2 3 2 5. 4. 4 0 0.5 0.5
```

Crystal

3.52387 3.52387 3.52387 90. 90. 90.

1 0.0 0.0 0.0

1 0.5 0.5 0.0

1 0.5 0.0 0.5

1 0.0 0.5 0.5

The configuration is $3d^9 4s^1$ with 5 electrons 3d up and 4 electrons 3d down, 0.5 electron 4s up and 0.5 electron 4s down.

In case of an antiferromagnetic structure, one must put a minus sign in front of the atom type number to specify the atom with the reverse spin. For example:

Atom

23 2 3 2 2. 0. 4 1 0.5 0.5

8 0

Crystal

7.255 5.002 5.548 90.0 96.75 90.0

1 0.34380 0.00080 0.29910

1 0.65620 0.99920 0.70090

-1 0.84380 0.50080 0.79910

-1 0.15620 0.49920 0.20090

2 0.40700 0.84500 0.65200

2 0.09300 0.84500 0.34800

2 0.59300 0.15500 0.34800

2 0.90700 0.15500 0.65200

Alternatively, it is possible to use the keyword « *Atom_conf* » which allows to specify the orbital occupancy of only some of the atoms:

Atom_conf

4 1 2 2 3 2 2. 0. 4 1 .5 .5 → nbr of atoms with the configuration, then index of these atoms

4 3 4 2 3 2 0. 2. 4 1 .5 .5 → then nbr of orb, then n,l pop up, pop down for each orbital...

Crystal

7.255 5.002 5.548 90.0 96.75 90.0 / a, b, c, alfa, beta, gamma

23 0.34380 0.00080 0.29910 V8

23 0.65620 0.99920 0.70090 V6

23 0.84380 0.50080 0.79910 V4

23 0.15620 0.49920 0.20090 V5

8 0.40700 0.84500 0.65200 O 1 8f

8 0.09300 0.84500 0.34800 O 2 8f

3) Spin axis

By default the spin axis is along the *c* axis but in the trigonal unit mesh in which it is along the *c* axis of the associated hexagonal unit mesh. This axis can be orientated along any (but uniform) direction:

Axe_spin

-0.08909 0. -0.15025 → in unit mesh

Another way to specify this axis is using the Euler angles. The orthogonal basis to do that is such that \mathbf{z} is along \mathbf{c} , but for trigonal system in which \mathbf{z} is along the hexagonal axis, \mathbf{x} is along $\mathbf{b} \times \mathbf{c}$ and \mathbf{y} is along $\mathbf{z} \times \mathbf{x}$. One then has to write:

Ang_spin

45. 90. 0. → rotation around \mathbf{z} , then around new \mathbf{y} , then around new \mathbf{z} . In this case spin axis is along (1,1,0) in the internal basis

Note that it is possible to define non collinear spin using atomic local basis. In this case one has to define the Euler local angles after the keyword “*crystal*” or “*molecule*”, in the line just before the corresponding atoms. See “non-spherical atom”. Example:

Crystal

7.7400 7.7400 3.8400 90. 90. 90.

45. 90. ! spin axis along (1 1 0)

1 0.0 0.0 0.0

2 0.25 0.25 0.5

135. 90. ! spin axis along (-1 1 0)

1 0.5 0.0 0.0

2 0.75 0.25 0.5

315. 90. ! spin axis along (1 -1 0)

1 0.0 0.5 0.0

2 0.25 0.75 0.5

225. 90. ! spin axis along (-1 -1 0)

1 0.5 0.5 0.0

2 0.75 0.75 0.5

4) Self-consistent calculations

By default a calculation cycle is performed over the occupied states to determine the Fermi energy. The XANES calculation corresponds to the cycle 2, but this one is done with the initial potential. To avoid this Fermi (non SCF) energy determination, use the keyword :

No_Fermi

To perform a self-consistent calculation, use the keyword:

SCF

The calculation is then really self-consistent. The potential is calculated again at each cycle. The final XANES calculation uses this SCF potential. By default the convolution uses the calculated Fermi energy as cutting energy.

When using SCF calculation, the cluster radius for the SCF is the same than the one for the XANES. For a single evaluation in one cycle of the Fermi energy the radius is equal to the

minimum between the XANES cluster radius and 3.5 Å. In both cases, this radius can be changed using the keyword :

R_self

3. → Value of the radius for the SCF and/or Fermi energy evaluation.

With SCF calculation, the maximum number of cycle is fixed to 30. To modify this number use the keyword:

N_self

20

The electronic density of the cycle $n+1$ is interpolated between the electronic density of the two previous cycles. The initial weight of the cycle is 0.1. This weight decreases or increases according to the quality of the convergence between 4 and 0.25 times the initial value. To change this initial weight use:

P_self

0.05 → Initial weight.

The convergence is considered as realized when the variation of the total energy is less than $N_a \times \Delta E_{\text{conv}}$, in which N_a is the number of atom and ΔE_{conv} a criteria fixed at 0.1 eV. To modify ΔE_{conv} put:

Delta_E_conv

2. → Value of ΔE_{conv} .

It can happen, that when the weight increase up to 2 or 4 times the p_{self} value, the calculation become very unstable (it can be seen on the orbital occupancy). One can then limit this increase by the parameter $p_{\text{self_max}}$:

P_self_max

0.05 → value of the maximum weight. When used, the best is to put it equal to p_{self} .

The SCF calculation is done with a non-excited cluster. The excited potential of the absorbing atom is taken into account only in the K and L_1 edges in the XANES final calculation by difference. If one wants to perform a SCF calculation with an excited cluster, use the keyword:

SCF_exc

If one wants perform a simulation with SCF in excited state but without any screening, that is with a cluster charged +1, with SCF_EXC one has to use keyword “*Screening*” with the value 0 (see this keyword).

For a magnetic calculation, by default the spin polarization is kept fixed in amplitude. The total number of spin up and spin down electron is fixed along the self-consistent (for an anti-ferro, for the total on the atoms, the number of majoritarian spin and minoritarian spin electron are kept fixed). To have it free (equivalent to version before 7th of June 2012), use the keyword:

SCF_mag_free

5) Hubbard correction

To include a Hubbard correction, use the keyword:

Hubbard

5. → value in eV of the Hubbard parameter (U-J).

This correction is applied to the “localized valence orbitals” of the chemical specie first in the list under “*crystal*” or “*molecule*” or under “*atom*” if used. Most often, it is convenient to perform also a self-consistent calculation. To apply different Hubbard corrections on different atoms, just put several numbers in the same line, following the order of the chemical species. Note that in principal this correction can be applied only on insulating compounds.

6) TDDFT calculations

TDDFT calculations can be performed using the keyword:

TDDFT

Then a first mono-electronic (LSDA) calculation is performed, followed by the TDDFT cycle. Then a correction is performed making a mixing between the edges. This can be useful for the L₂₃ edges of the transition 3d elements. Calculations use a local Kernel with its Coulomb and exchange-correlation parts. To have only the Coulomb part, use the keyword:

RPA

The TDDFT simulation contains a matrix inversion whom dimension is proportional to the number of spherical harmonics characterized by a l maximum equal to lmax. In optics lmax = 2 when Z < 21 and 3 in the other cases. To save time, especially in optics, it can be useful to limit it with the keyword:

Lmax_tddft

2 → maximum l value for this part of the calculation.

7) Conter-ion and Helmholtz layer

It is possible to add a potential simulating the effect of counter-ions in the case of a charged molecule imbedded in a solvent or to consider a non-ordered Helmholtz layer on top of an electrochemical interface (see the section on 2D diffraction). The corresponding potential is given by : $\Delta V = \frac{\sqrt{\pi}}{2} V_{helm} \operatorname{erf}\left(\frac{z-z_0}{\alpha \Delta_{helm}}\right) / \left(\frac{z-z_0}{\alpha \Delta_{helm}}\right)$ where z is the radial coordinate or the coordinate perpendicular to the surface, z₀ is the distance between counter_ion (or Helmholtz layer and most outer atom and Δ_{helm} is the half height width, with avec α = 0.285925223. When finite different method is used, this potential is also added on all the grid of point inside the interval. All this is obtained with:

Helmholtz

-10. 2. 8. $\rightarrow V_{helm}, \Delta_{helm}, P_{helm}$ (en eV et Å)

If, Δ_{helm} is omitted, its default value is : $\Delta_{helm} = 2z_0$, which corresponds to an added potential on the most outer atom equal to $\frac{1}{2}V_{helm}$. Alternatively one can use another potential shape for $z \in [z_0, z_0 + P_{helm}]$ with $\Delta V = \frac{1}{2}V_{helm} \left(1 + \cos\left(\pi \frac{z-z_0}{\Delta_{helm}}\right)\right)$ and $\Delta V = 0$ outside. For is write:

Helm_cos

-10. 2. 3. $\rightarrow V_{helm}, z_0, \Delta_{helm}$ (en eV et Å)

8) Exchange-correlation potential

By default the calculations are done using the real Hedin, Lundqvist and Von Barth potential. If one wants to use the Xalpha potential, one must introduce the keyword "*Xalpha*" followed by the value of the corresponding parameter:

Xalpha

0.3333 \rightarrow value of the Xalpha parameter

Remark : the value of the Xalpha parameter the closest to the other models is 2/3.

In the case of using the exchange-correlation potential coming from "*flapw*", if one wants to keep this potential independent of the energy put the keyword "*Xalpha*" with any value beneath. The potential will not be Xalpha but the one calculated by *flapw*.

If one prefers to use the Perdew and Wang's potential, put the keyword:

Perdew

This potential will be used with the parametrization of Moruzzi, Janak and Williams but keeping the energy dependency proposed by Hedin and Lundqvist.

9) Screening

By default, the screening is one electron on the first non-full valence orbital of the absorber. If one wants to modify this value or the number of the orbital, use the keyword "*Screening*" followed by the quantum numbers of the valence orbital and the value in number of electron (better inferior to 1 !) of the screening. In this case the cluster is not anymore necessary neutral.

Screening

3 2 0.2 \rightarrow n, l, and screening on the valence orbital

One can omit the orbital quantum number. In this case the code, by default will choose the first non-full valence orbital and one just write:

Screening

0.2 \rightarrow Screening on the valence orbital

For magnetic calculations, screening is by default the half of the total screening electron on both spin states. Note that the occupancy rate cannot exceed $2\ell+1$ for each spin state. If it is the case the exceeding amount is set on the other spin state. Whatever, one can impose a different screening, polarized in spin, just writing 2 numbers below « *Screening* », eventually following n and ℓ :

Screening

0.5 0.3

→ screening 0.5 up, 0.3 down

VI- Sophisticated keywords

1) Use of atomic electronic densities coming from external files

When using electronic densities coming from external files, these ones must be specified with the keyword “*atom*”. One eventually has to specify one atom more because the absorbing atom is a special chemical specie different of the non-excited atoms. Its files must also contain the electron initial orbital before excitation (for instance the 1s one for a K edge). For the format of these files see chapter F. In the case there are two atoms:

<i>Atom</i>	→ keyword preceding the atomic electronic densities
Fe_exc.txt	→ name of the file containing the excited atom
Fe.txt	→ name of the file containing the type 1 atom
Oxygen.txt	→ name of the file containing the type 2 atom

One can also put one or several atoms calculated internally among the exterior files:

<i>Atom</i>	
Fe_exc.txt	→ name of the file containing the excited atom
Fe.txt	→ name of the file containing the type 1 atom
8 2 2 0 2. 2 1 4.	→ type 2 atom calculated internally

Once defined the atoms types, the molecule or mesh description uses the atom type number and not the atomic number. For example, for the cfc copper crystal:

<i>Atom</i>	
c:\users\joly\dirac\cu\trdrslcu.e01	→ excited copper file
c:\users\joly\dirac\cu\trdrslcu.001	→ standard copper file

<i>Crystal</i>	→ Crystal structure
3.610 3.610 3.610 90. 90. 90.	→ a, b, c, α , β , γ
1 0.0 0.0 0.0	→ Type number, position
1 0.5 0.5 0.0	
1 0.5 0.0 0.5	
1 0.0 0.5 0.5	

The first index is now the atom type number (but the excited one) in the list under “atom”.

2) Uses of atomic electronic densities of Clementi and Roetti

By default the program calculate the atomic basis using a Hartree-Fock-Dirac procedure. If you prefer the Clementi and Roetti basis use the keyword:

Clementi

Note that the Clementi basis exists only for $Z < 55$. When using this basis, the self-consistent calculations are not anymore possible as well as the automatic evaluation of the Fermi level (The program takes automatically the option « *No_Fermi* »).

3) Non spherical atomic electronic densities

When calculating using the finite difference method, it is possible to define non-spherical atomic electronic densities. This is done with the keywords "*Atom_conf*" and "*Atom_nsph*". For each non-spherical atom one performs an expansion in spherical harmonics of these non-spherical orbitals. In practice, one writes an expansion in (l,m) of each of these orbitals with at the end the number of electron it contains. Optionally a rotation of the local atomic basis can be performed using the Euler angles (see keyword "*Ang_spin*", for the Euler angle definition). In this case, the Euler angles must be given in the line in front of the atom. Example:

Atom_conf

```
1 1 1 3 2 3.
```

Atom_nsph

```
1 1 1 ! nbr of atoms, index of the atoms, nbr of orbitals
0. 0. 0. 0. 0. 0. 1. 0. 0. 1. → (l,m) in natural order, here (l=2,m=0) has 1 electron
```

Molecule

```
1.9 1.9 1.9 90. 90. 90. → a, b, c, α, β, γ
0. 45. 0. → 45° rotation around the Oy axis
23 0.0 0.0 0.0
8 1.0 0.0 0.0
8 -1.0 0.0 0.0
8 0.0 1.0 0.0
8 0.0 -1.0 0.0
```

Note that the code does not take into account the eventual decrease in symmetry induced by this configuration. One thus must take care.

4) Mesh or molecule charge

The unit mesh must be neutral. A molecule is also often neutral. A test is performed in the program to verify this neutrality. If one wants to omit it, put:

Chlibre

5) Orbital dilatation

It is possible to modify the valence orbitals defined above expanding or contracting them. This can be very useful for ionic material, for instance oxygen 2^- , in which the atomic bases are calculated for neutral atoms. For this purpose introduce "*Dilatorb*" then for each orbital one wants to dilate, the atomic type (number in the list "*Atom*" or "*Atom_conf*"), the valence orbital number and the expansion coefficient. For the absorbing atom, add the same lines, but with "0" in first column, in order both excited and non-excited atoms are expanded:

Dilatorb

```
3 2 0.3
0 2 0.3    → for the "excited absorbing atom"
```

6) Reference of the photoelectron wave vector

By default the reference of the electronic wave vector (and so of the kinetic energy outside the sphere of calculation and the muffin-tin ground in case of muffin-tin calculation) is taken as the average of the potential between the absorbing atom and the first crown when calculating in green and at the outer sphere frontier in FDM. It is possible to impose this ground potential at the Fermi energy. The dependence versus the kinetic photoelectron energy is then added automatically to this term. Put then the keyword:

V0imp

```
-11.5    → value of the ground potential at the Fermi energy.
```

7) Maximum potential

For a calculation performed on a molecule, the potential increases when going away from the atoms. When this molecule is not in the gaz phase, but in solution or in any relatively dense surrounding, this increase is artificial. In case of calculation under the finite difference mode, this reach to false bounded states very thin in energy. It can be useful to give a maximum value to this potential to avoid this phenomena with the use of the keyword "*Vmax*", followed by the value of this potential :

Vmax

```
-6.
```

8) Complex energy

It is possible to use a complex energy when working in the multiple scattering mode (green). For this, one has to specify the imaginary part (positive) of this energy in a table versus the photoelectron energy, under the keyword "*Eimag*":

Eimag

```
0. 0.5
10. 0.7
30. 3.
```

50. 5.
100. 6.

When a uniform broadening is sufficient, it is not necessary to specify the photoelectron energy:

Eimag
0.1 → value of the uniform width (eV)

The use of a small width (0.1 eV) is sometimes useful for the calculations at low kinetic energy of the photoelectron because the localized level (3d or 4f) can be too thin in energy to be correctly evaluated. It is even truer for photoemission calculation (see keyword "*photoemission*"). If this broadening is small it has no effect on the forthcoming convolution to take into account the widths of the hole and final states.

9) Radius of the cluster for the superposition of the potential

The potential inside the sphere of calculation with a radius R_s set under the keyword "*radius*", is calculated by superposition. To avoid frontier problems, the atoms taken into account for the superposition are all the atoms inside the calculation sphere plus an outer shell. By default this outer shell is 2.5 Å thick, giving thus a new sphere with a $R_s + 2.5$ radius. If one wants a bigger radius (useful for the oxides) put the keyword:

Rpotmax
15. → radius of the cluster for superposition in Angstrom

10) $\ell+1$ approximation for the selection rule

If one wants to make the $\ell + 1$ approximation in which only the transition $\Delta \ell = +1$ is authorized (so one neglects $\Delta \ell = -1$, for example for the threshold L_2 and L_3 , the transitions are only toward the *d* states), put the keyword:

lplus1

In the same way if you want only the $\ell - 1$ states put:

lminus1

11) Rydberg series

One can add an outer sphere having a $-1/r$ potential to analyze Rydberg series. In this case put the keyword "*Rydberg*":

Rydberg

12) Cluster origin

By default, the cluster origin is set on the absorbing atom. If one wants that this origin is not changed from the indata file put the keyword:

Noncentre

If one wants to impose a specific center put:

Center

0 0.25 1 → coordinates in unit cell parameter unit of the center

The center can be automatically chosen at the center of the cluster when writing nothing after the keyword center:

Center

The center can be chosen at the center of the cluster considering only the chemical specie of the absorbing atoms with the keyword:

Center_abs

To have the center on an axis parallel to c , (usefull for surfaces) write the keyword:

Center_s

0.124 0.5 → coordinates (x,y) in unit cell parameter unit of the axis

13) Cut of the potential

To get a cut of the potential in the output, along a line or a plan inside the FDM grid of points (does not work when using the multiple scattering theory (Green)), put:

Trace

1 1. 0. 0. 0. 0. 0. → n, a, b, c, p, q, r
 n = 1, cut along a line with vector (a,b,c) crossing the point (p,q,r) in unit mesh parameter
 n = 2, cut along the plane $ax + by + cz = p$
 n = 3, gives all the points

14) Density of state

To get the partial projection of the density of states and their integral for any harmonics on the absorbing atom, put the keyword:

Density

Then one gets a new output file with the extension *_sd0.txt*. It contains, in column, first the integral of the total atomic electron density, then the density and its integral of each (ℓ, m)

followed by the sum over m , that is the density and its integral for each ℓ . For magnetic calculation, the expansion is split between the “up” and “down” components. By default, real harmonics are used because they are directly the familiar p_x , p_y , p_z , d_{xy} , d_{xz} ... states. The correspondence is the following:

(0,0)	(1,-1)	(1,0)	(1,1)	(2,-2)	(2,-1)	(2,-0)	(2,1)	(2,2)
s	p_y	p_z	p_x	d_{xy}	d_{yz}	d_{z^2}	d_{xz}	$d_{x^2-y^2}$

It can be interesting, for example to get the orbital magnetic moment, to prefer output with the usual complex spherical harmonics. For this, use the keyword:

Density_comp

General remark on the values one gets: Usually for a transition metal, at the end of the d band, the integral must be close to 10 (because there are 10 d electrons!). At the end of the rather unlocalized s and p bands, the integrals rarely reach 2 or 6, because the electrons are counted only inside the atomic sphere of radius R_{mstd} (see in the `bav` file for its value). These orbitals having a rather big radius they are not all inside such atomic sphere. Calculation being performed in the continuum one also has to recall that n is not anymore a “good” quantum number. Consequently the integral continue to increase indefinitely with energy. At the end of a d (or f) band, rather localized, one nevertheless reaches more or less at 10 (or 14). In this case one can find the Fermi level energy just looking the energy in which the integral reach the supposed number of electron in the corresponding level of the atom. Sometimes the integral never reach 10 (or 14) or goes far higher. One of the reasons can be that the starting energy is too high, one thus loses the beginning of the band, and another reason can be the energy grid is not sufficiently thin. One has, in this case to take a path smaller (down to 0.01 eV or even less for $4f$ elements). In order to avoid such a thin grid, it is also possible to broaden a bit these localized states using the keyword “*Eimag*” (this works only in “Green” mode).

To get the projection on all the atoms in the cluster (or unit cell) put the keyword “*Density_all*” in place of “*Density*”. Then one gets n new output files with the extensions `_sdi.txt`, “ i ” being the atom number:

Density_all

15) Spherical tensors

To get the spherical tensors (in number of electron) put the keyword:

Spherical

One gets new files with the extension `_sph_atom1.txt` and `_sph_xtal.txt` for the atom and crystal. If one uses such file as indata for the convolution process, with the keyword “*fprime*”, one gets the f and f' of each tensor component of the atom. Another file `_sph_int_atom1.txt` contains the integral of the spherical tensors.

To get the contributions of the tensor components on each polarization and reflections, put in place of “*spherical*” the keyword:

Sphere_all

Then one gets a series of output files with the extensions *_sph_signal_xan.txt*, *_sph_signal_pol1.txt*, *_sph_signal_rxs1.txt* for the contributions on each polarization and reflections. The number after *pol* or *rxs* is the number of the polarization or reflection.

16) Cartesian tensors

To get the atomic cartesian tensors put the keyword:

Cartesian

One gets a new output file with the extension *_car.txt*. Then, one can use other keywords making that from this output file will be calculated f'' and f'' for any tensor component. For this purpose put the keyword "*fprime*" in the indata file.

17) Calculation area boundary

By default in FDM, the meshing is performed in a sphere extending up to the last atom inside the sphere of radius given under "*Radius*" plus the atomic radius (by default 0.65 Å) plus one inter-point distance (0.2 Å by default). In order to use a bigger sphere use:

Overad

1.2 → distance over the last atom + its radius to take into account.

18) Displacement of the absorbing atom

To move the absorbing atom in reference to its position given under "*molecule*" or "*crystal*" use:

Dpos

0.2 0.0 0.0 → displacement vector in Angstrom

19) Getting back the tensor components from a previous calculation

It is possible to extract the tensor component from the "bav" file of a previous calculation. At this step, the tensors must had been written thanks the use of "check_coabs" keyword with the number 2 below. To get them one must write:

Extract_tens

Sim/calcul_prec_bav.txt

In case of calculations for different absorber position whose number is the number of output files, one has to extract the tensors from these n files if there are n nonequivalent absorbing atoms. This is automatic. If the calculation is to be performed on only one or some of the non-

equivalent atom sites, this must be specified by the keyword "*Absorber*" followed by the keyword "*extractpos*" to select in the bav file the corresponding atom:

Absorber

4 8

Extractpos

2 3 → in the bav file, the prototypical atoms 2 and 3 are chosen. Their tensors are used respectively by the atoms 4 and 8 in the list in the indata file.

It is possible to make rotations of the local repair. It can be useful when using tensor calculated in a different base. It is possible de make 1, 2 or 3 rotations around the axes Oz, then Ox, then Oy. Each rotation defines a new local base, starting point of the next rotation. For this write:

Rotsup

30. 0. 45. → 30° rotation around Oz, 0° around Ox then 45° around Oy.

For the same purpose on can apply symmetry operation on the extracted tensors. They are codified as explained in the manual with the "*symsite*" keyword. They must be given in the same order than the extracted files. These codes are given after the keyword *extractsym*:

Absorber

4 8

Extractsym

1 25 → code 1 for identity, 25 for inversion.

When keyword "*Extract*" or "*Extract_ten*" are used, no "bav" file is created, but when one also uses keyword "*check*".

20) Energy shift of the spectra

If one has gotten a reference for the initial orbital, it is possible to give it under the keyword "*Epsii*". This will produce a shift of the output spectra equal to the difference between this energy and the energy calculated in the program. It is safer to perform this operation with the shift parameters during the convolution step.

Epsii

6253.1 → positive value in eV.

21) Spectra by atom

If one wants to get the signal coming from each atom one has to add the keyword:

Allsite

One then gets, at most of the usual output files, new output files of the type atom1.txt, atom2.txt, etc... The suffixes *_atom1*, *_atom2* correspond to the number of the atom.

22) Use of densities and potential coming from FLAPW

If the potential is imported from a FLAPW calculation, one has to introduce the keyword "*flapw*" in place of "*atom*". The keyword "*crystal*" or "*molecule*" becomes also unnecessary. Then must stand the names of the 5 FLAPW output files in case of spin un-polarized calculation and the 7 files in case of spin polarized calculation:

<i>Flapw</i>	→ names of the output FLAPW files
<i>tio2.struct</i>	structure and symmetry
<i>tio2.vcoul</i>	coulombian potential
<i>tio2.r2v</i>	exchange-correlation potential
<i>tio2.clmsum</i>	electronic density
<i>tio2.tils</i>	initial wave function

In case of polarized calculation the keyword "*magnetism*" must be before the keyword "*flapw*":

Magnetism

<i>Flapw</i>	
<i>tio2.struct</i>	structure and symmetry
<i>tio2.vcoul</i>	coulombian potential
<i>tio2.r2v</i>	exchange-correlation potential spin up and spin down
<i>tio2.clmsum</i>	total electronic density
<i>tio2.clmup</i>	valence electron electronic density, spin up
<i>tio2.clmdn</i>	valence electron electronic density, spin down
<i>tio2.tils</i>	initial wave function

The last file contains the wave function of the initial core orbital. If one is ok with the initial wave function calculated internally, one can avoid it. One must in this case substitute the keyword "*flapw*" by "*flapw_psi*".

If one wants to have an energy dependent exchange-correlation potential, one must add the keyword "*Hedin*".

By default, the absorbing atom is the first one in the "*struct*" file list. If one wants that it is the n^{th} , put the keyword:

absorber
n

For the actual version of wien2k and magnetic calculations, the exchange-correlation potential for spin up and spin down is given in two different files. In this case one must use the keyword:

<i>Flapw_n</i>	
<i>tio2.struct</i>	
<i>tio2.vcoul</i>	
<i>tio2_up.r2v</i>	exchange-correlation spin up
<i>tio2_dn.r2v</i>	exchange-correlation spin down
<i>tio2.clmsum</i>	

tio2.clmup
tio2.clmdn
tio2.tils

when not using the wien2k core state use :

Flapw_n_p
tio2.struct
tio2.vcoul
tio2_up.r2v
tio2_dn.r2v
tio2.clmsum
tio2.clmup
tio2.clmdn

Note also that sometimes, it is necessary to increase the value of number of point in the radial mesh. In this case the value to put is given in an error message. It is possible using the keyword *nrato*.

Nrato

950 → new value given in the error message

23) Jump in the indata file

It is possible to jump over a part of the information written in the indata file with the keywords "*Jump*" and "*Endjump*". All what is between them is not read. For example:

Jump

Quadrupole
Spinorbit

Endjump → the keywords "*Quadrupole*" and "*Spinorbit*" are not considered.

24) Comment

It is possible to introduce a line of comment which will be copied in the output files. For this purpose write:

Comment

Iron K-edge in magnetite

→ line of comment

25) Atomic spectra

To get in last column the atomic absorption spectra (without the neighbour atoms), put the keyword:

Xan_atom

26) Different absorbing atoms calculated in one run

The electronic structure is calculated in all the cluster of calculation and thus in all atoms in it. Consequently, it is in principal possible to get the absorption spectra of all the atoms in only one run. In principal, the absorbing atom is nevertheless “excited”, thus becomes intrinsically different and one calculation must be performed for each absorbing atom. When neglecting this difference, it is possible to get all the absorption spectra of the all the atoms, equivalent and non-equivalent by symmetry operation, of the same chemical specie in only one run using the keyword:

One_run

For the L₂₃ and M₄₅ edges, the default calculation is such that the absorbing atom is not excited, thus the result in *one_run* is in principal good. With “*One_run*” all the absorbing atoms are calculated, one thus must take care that they are all in the sphere of calculation. The use of keyword “*Center*” or “*Center_abs*” is sometimes helpful in this context.

27) Calculation of the Crystal Overlap Occupation Population (bounding)

The overlap between orbitals of neighboring atoms can be obtained using the keyword:

COOP_atom

1 2 4 → indices of the selected atoms

The selected atoms are the ones around which one looks for the first neighbors up to a default distance equal to the sum of the atomic radii. It is between these couples of atom that are calculated the orbital overlaps.

These ones are given by the formula $n_{COOP} = S_{ab} \int \varphi_a(\mathbf{r})\varphi_b(\mathbf{r})d\mathbf{r}$, where S_{ab} is the overlap integral and a and b , 2 neighboring atoms. n_{COOP} , positive or negative gives the bounding and antibonding. Results are in files of the type « *_coop_n_m* » where n and m are the indices of the 2 neighboring atoms. The files contain versus energy, n_{COOP} total, followed by $n_{COOP}(\ell_a, \ell_b)$, then $n_{COOP}(\ell_a, m_a, \ell_b, m_b)$.

We can modify the limit distance between 2 neighboring atoms to search for second neighbors for instance. This is done with:

COOP_dist

4 → Maximum distance in Angström

VI- Technical keywords

Here stand the technical keywords necessary for specific tests and optimization of the code, for expert users.

1) Details on the calculations

To get details on the calculations done in every routine, put the keyword:

Check

3 3 3 3 1 1 1 2 1 1
1 1 2

Up to 30 numbers can be given. The correspondence between the number and the routine is given by :

1 : Lecture	2 : Atom / Dirac	3 : Symsite	4 : Init-run
5 : Agregat	6 : Polarization	7 : Atom_selec	8: Etafin
9 : Reseau	10 : Potato	11 : Orbval	12 : Pot0
13 : Potentiel	14: Enrgseuil	15 : Ylm	16 : Potex
17 :	18 : Sphere	19 : Mat / MSM	20 : Tenseur
21 : Coabs	22 : Tddft – sphere	23 : Tddft – Kernel	24 : Tddft - Chi_0
25 : Tddft – Chi	26 : Hubbard	27 : SCF	28: State
30 : Convolution			

Values can go from 0 to 4 giving more and more for each routine. By default, there is 1 for all the routines. One can also write:

check_all	→ equivalent to icheck = 3 for all the subroutines
no_check	→ equivalent to icheck = 0 for all the subroutines
check_conv	→ equivalent to icheck = 3 for the convolution
check_pot	→ equivalent to icheck = 3 for all the subroutines concerning the potential
check_mat	→ equivalent to icheck = 3 for the matrices MSM or FDM
check_sph	→ equivalent to icheck = 3 for the subroutine “sphere”
check_coabs	→ equivalent to icheck = 3 for the subroutine “coabs”
check_tens	→ equivalent to icheck = 3 for the subroutine “tenseur”

2) Symmetry

The point symmetry is calculated automatically. Anyway the expert users can impose it by the keyword:

Sym

2/m → Schoenflies coefficient or international table of crystallography coefficient. Be careful that these symmetries are possible only for the mesh axis defined in the indata.

3) Calculation basis

By default the connection to the continuum is performed in real base (neuman, bessell). To use a complex base (bessel, hankel) put:

Basecomp

4) Equivalent atoms

The program automatically calculates the total signal resulting from the equivalent atoms by symmetry. It is nevertheless possible, for the expert *fdmnes* user, to impose the symmetry relation between the atoms and thus to impose this summation. One must then use the keyword "*symsite*". Following this keyword, stands the number of non-equivalent atoms (or number of group of atoms), then for each of them, the number of equivalent atoms followed the list of the relative symmetry to the first atom of the list and the atomic position in mesh unit. The symmetry are codified by number going from 1 to 64, see the list given bellow.

Symsite

```

3          ! Number of non-equivalent atoms (or group of atoms)
4          ! Number of equivalent atoms, group 1
 1 0.2500 0.2500 0.2500
24 0.7500 0.7500 0.2500
23 0.7500 0.2500 0.7500
22 0.2500 0.7500 0.7500
4          ! Number of equivalent atoms, group 2
 1 0.0000 0.0000 0.0000
22 0.0000 0.5000 0.5000
23 0.5000 0.0000 0.5000
24 0.5000 0.5000 0.0000
4          ! Number of equivalent atoms, group 3
 1 0.5000 0.5000 0.5000
23 0.0000 0.5000 0.0000
22 0.5000 0.0000 0.0000
24 0.0000 0.0000 0.5000

```

The negative indices correspond to the same symmetry plus time reversal (in case of magnetism).

Symmetry code:

1: identity	7: rot $4\pi/3$ around (-1,1,1)
2: rot $2\pi/3$ around (1,1,1)	8: rot $2\pi/3$ around (1,1,-1)
3: rot $4\pi/3$ around (1,1,1)	9: rot $4\pi/3$ around (1,1,-1)
4: rot $2\pi/3$ around (1,-1,1)	10: rot $2\pi/2$ around (1,1,0)
5: rot $4\pi/3$ around (1,-1,1)	11: rot $2\pi/2$ around (-1,1,0)
6: rot $2\pi/3$ around (-1,1,1)	12: rot $2\pi/2$ around (1,0,1)

- | | |
|--|---|
| 13: rot $2\pi/2$ around (-1,0,1) | 38: rot $2\pi/3$ around (1,1,-1) and inversion |
| 14: rot $2\pi/2$ around (0,1,1) | 39: rot $4\pi/3$ around (1,1,-1) and inversion |
| 15: rot $2\pi/2$ around (0,-1,1) | 40: plane perpendicular with 0x |
| | 41: plane perpendicular with 0y |
| 16: rot $2\pi/4$ around 0x | 42: plane perpendicular with 0z |
| 17: rot $2\pi/4$ around 0y | |
| 18: rot $2\pi/4$ around 0z | 43: diagonal plane $y = z$ containing 0x |
| 19: rot $-2\pi/4$ around 0x | 44: diagonal plane $x = z$ containing 0y |
| 20: rot $-2\pi/4$ around 0y | 45: diagonal plane $x = y$ containing 0z |
| 21: rot $-2\pi/4$ around 0z | 46: diagonal plane $y = -z$ containing 0x |
| | 47: diagonal plane $x = -z$ containing 0y |
| 22: rot $2\pi/2$ around 0x | 48: diagonal plane $x = -y$ containing 0z |
| 23: rot $2\pi/2$ around 0y | |
| 24: rot $2\pi/2$ around 0z | 49: rot $2\pi/3$ around 0z |
| | 50: rot $4\pi/3$ around 0z |
| 25: inversion | 51: rot $2\pi/6$ around 0z |
| | 52: rot $10\pi/6$ around 0z |
| 26: rot $2\pi/4$ around 0x and inversion | 53: rot $2\pi/3$ around 0z, negative axe |
| 27: rot $2\pi/4$ around 0y and inversion | 54: rot $4\pi/3$ around 0z, negative axe |
| 28: rot $2\pi/4$ around 0z and inversion | 55: rot $2\pi/6$ around 0z, negative axe |
| 29: rot $-2\pi/4$ around 0x and inversion | 56: rot $10\pi/6$ around 0z, negative axe |
| 30: rot $-2\pi/4$ around 0y and inversion | |
| 31: rot $-2\pi/4$ around 0z and inversion | 57: plane at 30° containing 0z |
| | 58: rot $2\pi/2$ around axe at 30° perp. 0z |
| 32: rot $2\pi/3$ around (1,1,1) and inversion | 59: plane at 60° containing 0z |
| 33: rot $4\pi/3$ around (1,1,1) and inversion | 60: rot $2\pi/2$ around axe at 60° perp. 0z |
| 34: rot $2\pi/3$ around (1,-1,1) and inversion | 61: plane at 120° containing 0z |
| 35: rot $4\pi/3$ around (1,-1,1) and inversion | 62: rot $2\pi/2$ around axe at 120° perp. 0z |
| 36: rot $2\pi/3$ around (-1,1,1) and inversion | 63: plane at 150° containing 0z |
| 37: rot $4\pi/3$ around (-1,1,1) and inversion | 64: rot $2\pi/2$ around axe at 150° perp. 0z |

5) Bounded states

If one wants to make positive the kinetic energy for the calculation of bounded states (beneath the average potential), put the keyword "*E_out_min*" followed by the kinetic energy in eV:

E_out_min
0.05 0.5

The first value is for the atomic scattering, the second one only used in FDM is for the cluster scattering. Default values are respectively 0.2 and 1 eV.

6) Inversion of the matrix in the multiple scattering mode

By default one inverts the matrix $(\kappa^{-1} + i + G)$ in which κ corresponds to the Neuman-Bessel normalization. To avoid eventual convergence problems (not seen yet), one can prefer to calculate $(1 + G\tau)^{-1}\tau$. For this use the keyword:

Normaltau

7) Muffin-tin potential

If one wants to use the muffin-tin approximation in FDM put the keyword:

Muffintin

8) Non excited absorbing atom

At the K, L1, M1 edges the absorbing atoms is calculated excited that is with a hole in the core level and an extra electron in the first non-occupied valence state. To perform a calculation with a non-excited absorbing atom, that is with an absorbing atom without a hole in its core level and without screening put:

Nonexc

At the other edges, the absorbing atom is not excited. To make it excited put the keyword:

Excited

9) Modification of the grid of point parameters

By default the finite difference calculation is performed at order 4 with an inter point distance equal to 0.25 Å, the radius of the spherically symmetric area is around 0.65 Å (less for the light elements). To modify them write:

Rmt → to modify the muffin-tin radius
0.65

Iord → To modify the order of the Taylor expansion.
2

Adimp → to modify the interpoint distance
0.20

It is possible to use inter-point distance, δ , decreasing with energy. This can be useful in EXAFS, because at high photoelectron kinetic energy, lower δ are necessary to get the same precision. For this, write :

Adimp

0.24 100. 0.20 250. 0.16 400. 0.12 500. 0.08

→ $\delta = 0.24$ up to 100 eV, $\delta = 0.2$ between 100 and 250 eV,... $\delta = 0.08$ beyond 500 eV

Remark: when using different cluster radius, the energy values in which δ change, can be the same than for the radius (see « *Radius* » keyword).

10) Expansion in spherical harmonics

An expansion in spherical harmonics is performed in the atoms both in the multiple scattering mode and in the finite difference mode, but with a smaller radius. The maximum value of ℓ is obtained from the formula $kR = \sqrt{\ell_{max}(\ell_{max} + 1)}$, in which k is the photo electron wave vector and R the muffin-tin radius (R_{mtg} in the *bav* file) in atomic units. By default, to this value one adds 1. To fix a value of ℓ_{max} independent of energy just write:

Lmax → The ℓ_{max} value is fixed at 3
3

If one wants to keep the energy dependent formula but adding 1, 2 or more to this formula, put a negative sign in front of 1 :

Lmax → The ℓ_{max} value is now given by the formula plus 2.
-2

When the energy is increasing, ℓ_{max} increases. By default the maximum value of ℓ_{max} is set at 2 for $Z = 1$, at 3 for $Z = 2$ then at 4 for $Z \leq 18$, then at 5 for $Z \leq 36$, then at 6 for $Z \leq 54$, then at 7 for $Z \leq 86$ and 8 for Z over. To avoid this limitation put the keyword:

Lmaxfree

When working using the finite difference mode, the connection to the outer sphere also needs an expansion in spherical waves. This one uses the same formula that for the atom but R is now the radius of the outer sphere. To modify this ℓ_{max} one proceeds exactly in the same way with positive value to fix an energy independent value, and negative value to keep the formula but adding to the formula. By default the additional value is 5. The key word is now *lmaxso*:

Lmaxso
15

It is possible to give a maximum value to the ℓ_{max} of the outer sphere writing:

Lmaxso_max

28 → maximum value

11) Muffin-tin radius

By default the muffin-tin radius (used in multiple scattering theory, for the electronic density calculation and for the SCF) is calculated to have a jump of potential the smallest as possible and closed between the atoms. One can modify it by the use of 3 different keywords.

If one wants that this radius is calculated using the Norman procedure introduce the keyword:

Norman

If one wants that the muffin-tin radius is the half value of the interatomic distance, introduce the keyword:

Raydem

By default there is a 10% overlap between the spheres. If one wants to modify this put the keyword "*overlap*" followed by the value of the overlap:

Overlap

0.15 → 15% overlap

If one wants that the muffin-tin radius be determined by the value it gets when the atomic potential is equal to the interstitial potential, put the keyword "*Rmtv0*", followed by the value of this potential (takes place the one defined by the keyword "*V0imp*"):

Rmtv0

-12.

If one wants to impose the radius for each chemical species put "*Rmtg_Z*" followed by the atomic number and the radius in Angstrom for each of them, and line after line:

Rmtg_Z

8 0.65 → Z = 8, Oxygen, radius = 0.65 Å

26 1.1 → Z = 26, Iron, radius = 1.1 Å

12) Continuity of the potential at the muffin-tin radius

To limit potential jump between the interstitial area and the muffin-tin sphere, a linear interpolation on the potential is performed in the last 10% of the radius. For this put:

Supermuf

13) Atomic charge

The program gives in the output file the atomic charge integrated up to the muffin-tin radius. If one wants to get this value for a different radius put:

Rchimp → corresponding keyword
1.2 → radius value for each type of atom
0.7

14) Line length

When reading a file, the maximum number of columns of number is by default 10001. This number can be insufficient for example when making resonant diffraction. It can thus been increased using the keyword:

Length_line
 20000 → maximum number of columns in one line

15) Resolution in initial state

The outputs are resolved in initial states only for magnetic calculation. For a non-magnetic calculation, it is possible to get anyway the output resolved be initial state using the keyword:

Core_resolved

It is then possible in the step convolution to select one of these states.

16) Memory save

For low symmetrical material, it can be useful to save some memory space. This makes an approximation (but a good one) on the potential calculation. For this use the following keyword:

Memory_save

C- Convolution

I- Introduction

The *fdmnes* program allows performing:

1) the convolution by a lorentzian of absorption spectra eliminating the occupied states. For the anomalous diffraction, it makes the integration over energy of the unoccupied states then calculates the intensity of the diffracted peaks.

2) a weighted summation over different outputs of the *fdmnes* program. This summation can include a relative shift between the spectra.

This step of the calculation can be done together with the previous one, just adding the corresponding keywords. It can also be done independently, writing another input file contains only the keywords related to the convolution. Here comes an example:

Calculation	→ To give the file name resulting from the previous step.
<i>g_rs43.txt</i>	→ name of the file
Conv_out	→ To give the file name of the convoluted spectra
<i>g_rs43_conv.txt</i>	→ name of the file
Convolution	→ keyword to specify an “arctangent” shape for the broadening”
E_cut	→ keyword to specify a cutting energy
<i>-1.5</i>	

II- Keywords for the convolution

1) File names

The indata files for the convolution step are the output files of the previous step whose names are defined under "*filout*". When the convolution is performed together with the main calculation, it is not necessary to specify again this name. In the other cases one has to introduce it or them with the keyword "*Calculation*":

Calculation	→ To give the file names resulting from the previous step.
<i>g_rs43.txt</i>	→ name of the file

Another example with 2 indata files and with different shifts:

Calculation	→ To give the file names resulting from the previous step.
<i>Fe_rs64_01.txt</i>	→ name of the file 1
<i>1.0 0.2</i>	→ weight, shift
<i>Fe_rs64_02.txt</i>	→ name of the file 2
<i>1.0 -0.2</i>	→ weight, shift

The weight is applied on the XANES and the anomalous scattered amplitudes. For the *DAFS* calculations, the output files, before convolution contain in second line the Thomson factors plus the anomalous part of the non-resonant atoms. By this way it is possible to take into account an average configuration for substitutions of atoms of different chemical specie on the same sites.

The name of the convolution output file is by default the input file name with the added suffix "*_conv.txt*". Anyway it is possible to impose another name with the keyword:

Conv_out

Fe_rs64_sum_conv.txt → name of the convoluted spectra file

To specify a working directory put the keyword *directory* followed by the directory name with at the end the separator (/ or \):

Directory

C:/Documents and Settings/joly/Mes documents/xanes/xanout/v2o3/

When there are more than one indata files, it is possible to have not only the total convoluted file but also the individual convolutions for all the indata files. For this write the keyword:

All_conv

Before the edge, the absorption is zero. It is possible to take into account the background coming from the edges of lower energy from all the chemical elements in the material. For this use the keyword:

Abs_before

2) Fermi level

The zero of the energy is set at the Fermi level. The states beneath the Fermi level are occupied. At the beginning of the convolution step, the cross section is thus set to zero, for negative energies. It is possible to change this cutting energy with the keyword:

E_cut → keyword to specify the new cutting energy (EFermi keyword in previous versions)
-1.5

The cutting value is applied without additional shift because in principle this last is used to simulate the core level shift. Anyway it is possible to apply it added with the *E_cut* value by the keyword:

Dec

It is also possible to have different values of the cutting energy for the different files. This is written in third column after each file name:

Calculation

Fe_rs64_01.txt

1.0 0.2 -2.

Fe_rs64_02.txt

1.0 -0.2 -1.9

→ To give the file names resulting from the previous step.

→ name of the file 1

→ weight, shift, cutting energy

→ name of the file 2

→ weight, shift, cutting energy

The cutting energy written under "*E_{cut}*" when specified, is then not considered.

3) Convolution width

The convolution to apply depends on the core level width (Γ_{Hole}) and the final state width. Γ_{Hole} is automatically determined. Nevertheless, it is possible to modify it with the keyword:

Gamma_hole

2.05


→ Core level width (Γ_{Hole}) in eV.

To have different values of Γ_{Hole} for a double edge like L₂₃, just put 2 numbers.

There are different ways to calculate the energy dependent broadening. The first one uses an arctangent formula:

$$\Gamma = \Gamma_{\text{Hole}} + \Gamma_m \left(\frac{1}{2} + \frac{1}{\pi} \arctan \left(\frac{\pi}{3} \frac{\Gamma_m}{E_{\text{Larg}}} \left(e - \frac{1}{e^2} \right) \right) \right)$$

with: $e = \frac{E - E_F}{E_{\text{cent}}}$ in which Γ_m , E_{cent} and E_{Larg} are respectively the maximum width (at high energy) of the final state, the center and the width of the arctangent function. The depth at the center of the arctangent is Γ_m/E_{Larg} . A typical curve is given below:

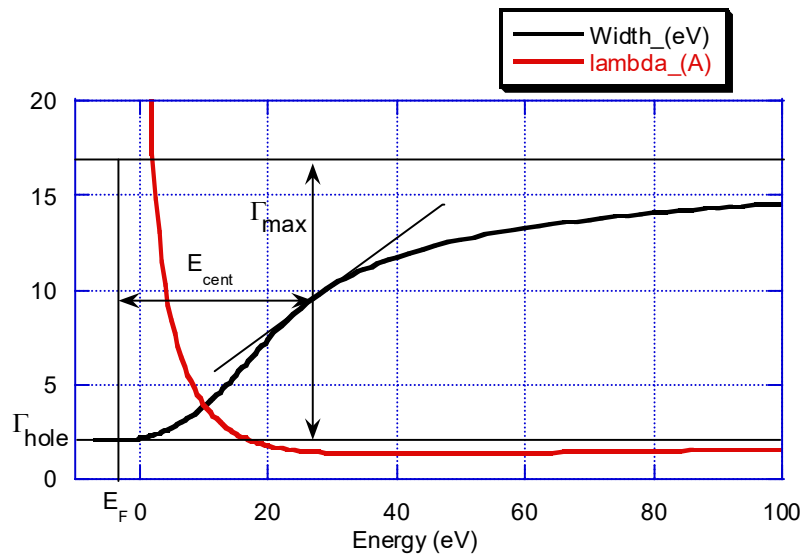

Such a convolution is introduced with the keyword:

Convolution

Default values for E_{larg} , E_{cent} and Γ_m are respectively: 30, 30 and 15 eV. It is possible to modify them with the keywords:

Ecent

30. $\rightarrow E_{\text{cent}}$

Elarg

30. $\rightarrow E_{\text{larg}}$

Gamma_max

20. $\rightarrow \Gamma_m$

In the convolution, along the integration it is the width of the running energy which is taken. It is possible to use the width of the final state energy corresponding to the energy of the elastic photon. One then makes the integration with a constant width. This procedure improves the agreement with experiment especially in the pre-edge range in which the other procedure increases the background. To impose nevertheless a variable width along the integration in XANES uses the keyword :

Gamma_var

It is also possible to use the Seah-Dench formula for the calculation of the broadening. In this case one gets:

$$\Gamma = \frac{A\Gamma_m E_p}{\Gamma_m + AE_p} + \Gamma_{\text{Hole}}, \quad \lambda = \frac{1}{A\sqrt{E_p}} + \frac{\sqrt{E_p}}{\Gamma_m}, \quad \text{with: } E_p = E - E_F.$$

This is performed with the keyword:

Seah

1. 20. $\rightarrow A, \Gamma_m$

One can also use a simple table with the keyword "*table*":

Table

0.0 1. \rightarrow Energy ($E-E_F$) and broadening $\Gamma(E-E_F)$
 10.0 1.5
 20.0 2.
 30.0 3.5
 100.0 5.

When the keyword "*E_cut*" is used without the keywords "*Convolution*" or "*Seah*", the width of the convolution is constant and equal to the width of the core hole. In all cases, the curve shape is written in the talkative "bav" file when the convolution is done in the same run than the main calculation. When convolution is done alone, to get this curve write the keyword:

Check_conv

To simulate the experimental resolution, it is also possible to convolute by a Gaussian. This convolution is performed after the lorentzian. For this write:

Gaussian → keyword
 1. 0.05 → width of the gaussian

This same convolution can be used as a first approximation to simulate the thermal vibration effects. This is valid only in the extended part of the spectra. In this case one applies the rule “ $kr = \text{constant}$ ” and the convolution width becomes proportional to the energy. The parameter to write is the relative vibration amplitude in relation to the interatomic distance. This is done, writing a second number under "Gaussian":

Gaussian → keyword
 1. → width of the Gaussian, relative amplitude of the vibration

4) Starting energy

To get output spectra starting at a lower energy put:

Estart
 -8. → Value of the starting energy (eV)

5) Photoemission

With the keyword "*photoemission*", it is the photoemission spectrum which is calculated. In such calculation the cutting is simply above the Fermi level and not beneath it as in xanes. The convolution broadening is the independent of energy and is fixed by default to the hole width. This with can nevertheless been modified by the use of the keyword "*Gamma_hole*". Note that would have interest having previously performed a calculation in "green" mode with a minimum energy width (keyword "*eimag*") at least 0.1 eV, all this within an energy range convenient for photoemission:

Photoemission

6) Damping due to multi-electronic effects

In a simple approach, the multi-electronic phenomena can be seen as a simple damping in the absorption cross section and in the anomalous scattering amplitude in diffraction. This damping not taken into account, by default, can be introduced using the keyword "*S0_2*". It is typically around 0.8.

S0_2
 0.8 → value of the damping

When setting $S0_2 = 0$, in diffraction, one keeps only the non-resonant term.

7) Selection of the transition from a peculiar initial state

The fdmnes output files, before convolution, contain the spectra by transition from all initial states when the calculations are magnetic and by threshold in the other case. It is possible to perform the convolution for only one, or some, of these components (avoiding thus their total summation) with the keyword "*Selec_core*" with below the selected state. The order of these states goes from $j_z = -j$ to j ; for example at a L_3 edge: -3/2, -1/2, 1/2, 3/2. For a L_{23} edge, L_2 is before L_3 . In the example below, applied to a K edge, the 1s state with $j_z = -1/2$ is selected. It is purely spin down and after transition probes the valence down states:

Selec_core

When the « *Mat_polar* » keyword in the main calculation has been used, in order to calculate the transmission through a sample, the sample thickness in μm^{-1} must be given. This is done with the keyword:

Sample_thickness

100. → Sample thickness in μm^{-1}

See also the "*Stokes*" and "*Stokes_name*" keyword, below, which also apply in absorption when the « *Mat_polar* » keyword have been used.

8) Thermal disorder

It is possible to take into account the thermic disorder using the Debye model. For this purpose use the keyword:

Abs_B_iso

0.72 → Debye parameter ($8\pi^2\langle u^2 \rangle$) in \AA^2

For the same purpose but using as Debye parameter $\langle u^2 \rangle$ instead of $8\pi^2\langle u^2 \rangle$ use:

Abs_U_iso

0.01 → Debye parameter ($\langle u^2 \rangle$) in \AA^2

When these parameter are given in the indata file of the main calculation (using "*Atom_U_iso*" or "*Atom_B_iso*"), they are replaced for this part of the calculation by these new ones.

III- Keywords specific for resonant diffraction

1) Azimuthal scan

If there is an azimuthal scan in anomalous diffraction, one has also to give the corresponding indata file names (type *_scan.txt*) and the new output file. This is set after the keyword "*scan_file*":

Scan_file

Fe_rs64_01_scan.txt → name of the indata dafs file 1

Fe_rs64_02_scan.txt → name of the indata dafs file 2

Scan_conv

Fe_rs64_scan_conv.txt → name of the new output file with the convoluted scan

When the convolution file is included in the "main" calculation file, it is not necessary to specify the scan file names because they are known by the code. In this case, if one wants to convolute also these files just uses the "*Scan*" keyword:

Scan

2) Thomson

Thomson term is calculated by default. It is possible to modify it using the following keywords.

To eliminate the Thomson factor f_0 (as for a forbidden beam) put:

Forbidden

It can be useful to impose specific value for each *DAFS* reflections for the non-resonant (Thomson) structure factors and the resonant one of the other atoms. That is: $\sum_a p_a (f_{0a} + f'_{0a} + if''_{0a}) \exp(i\mathbf{Q} \cdot \mathbf{R}_a)$, with f'_{0a} and f''_{0a} considered only for the other chemical species and in which p_a is the site a occupancy rate. This is the case for example for occupancy rate non equal to one. These complex terms taken as a constant in all the energy range can be introduced through the keyword:

Thomson

1.1670478E+02 1.0583769E+02 1.1477827E+02 -8.2670689E+01

In the second line stand these terms, here for 2 reflections. They replace the values calculated by the program and placed in the second line of the output file.

3) Getting f' and f''

To get in the output file, the f' and f'' values, put the keyword:

Fprime

In the output file for each reflection, there are three columns, f' , f'' then intensity. The intensity is for the mesh, f' and f'' are for one atom of the mesh. f' contains also f_0 . To get only f' , put also the keyword "*forbidden*".

It is also possible to get the atomic (but in the solid) f' and f'' values without *DAFS* simulations. For this the best is to run the main simulation with the keyword "*spherical*". In a second step one convolute with the keyword *fprime* the corresponding output file "..._sph_atom1.txt"

To get the isolated atomic f' and f'' :

Fprime_atom

The values of the atomic f' and f'' are then in the talkative "bav" file, eventually created at this time when the convolution is done in an independent job.

4) Corrections

It is possible to make a double absorption correction, corresponding to a measurement done at the maximum of the « rocking curve » and thus not integrated over it. For this write the keyword:

Double_cor

This correction is effective only with *self_abs* or *Full_self_abs* in the main calculation step. One gets then in the output file after convolution a new column with 'Id' which contains the intensity doubly corrected.

Another correction is possible to take into account a dead layer, that is a disordered layer at the surface, not participating to the diffraction, but damping the signal in the incoming and in the outgoing. This is done with the keyword:

Dead_layer

150 → thickness in Angstrom of the dead layer.

This correction is also done during the convolution part of the calculation.

The absorption correction depends on the incidence and exit angles with the surface. By default one uses the (rather good) approximation of taking the average angles. To use the exact angles, one must define the surface plane. This is done with :

Surface_plane

1 2 0 → indices of the surface plane

III- Getting other polarization conditions

One can be interested by material with birefringence or rotating power properties with effect in resonant diffraction or in absorption measured in transmission mode and with an analyser not necessarily parallel to the incoming polarization. One can make this choice on the analyser condition at the convolution step calculation when in the previous step, the main calculation had been performed with all $\sigma\sigma$, $\sigma\pi$, $\pi\sigma$ and $\pi\pi$ polarization conditions. In absorption calculation this is done when using the keyword “*Mat_polar*”.

In that case, to get the intensities (only for diffraction) corresponding to incoming polarization right and left, add the keyword:

Circular

It is also possible to get the intensity corresponding to any Stokes parameters ($S1 = +/-1$ sigma and pi, $S2$ diagonal $S3 = +1(-1)$ circular right (left)) or with or without analyzer, this one with any angle. For example without analyzer one writes:

Stokes

```
0 0 1      → S1, S2, S3
0 0 -1
0 0.31 0.95
0 0.31 -0.95
```

Stokes

```
0 0 1 30. 50.      → S1, S2, S3, rotation and Bragg analyzer angles
```

If the Bragg angle analyzer is not specify, it is supposed perfect, that is 45° :

Stokes

```
0 0 1 30      → S1, S2, S3, rotation analyzer angle (0 =  $\sigma$ , 90 =  $\pi$ )
```

It is possible then to specify column name in the output file to make simplest the column identification, writing:

Stokes_name

```
no_name      → default name
no_name
+00+31+95    → name appearing in the first line, in the corresponding column
+00+31-95
```

In absorption one can check the different contributions on the transmitted signal coming from birefringence or from dichroism alone. For this purpose, write:

Check_biref

Then for each polarization condition one gets 3 new columns in the output file with in their name “nB” for “no birefringence”, “nD” for “no dichroism” and “nO” for no off-diagonal components.

IV- Independent convolution by a gaussian

Independently from any calculation the fdmnes code makes possible to perform a simple convolution by a gaussian of a table of data. This table must have in first column the abscissa then a set of column with the ordinates. All the columns must have the same number of elements. These data can be of any type, for example experimental spectra. Any number of text lines can be before the data.

The indata file asking for such a convolution is of the type:

File_in → keyword demanding the file name containing the data to convolute

File_name_data → name of the file

Conv_gaus → keyword demanding the width of the gaussian

2.5 → width in the same unit that the abscissa in the first column of the data file

File_out → keyword demanding the output file name of the convoluted data

File_name_out_conv → name of the file

End

D- Parameter optimization

It is possible to compare the calculated spectra to the experimental ones with the help of metric distances and R factor. The metric distances to compare the individual theoretical ($f_{th}^{(i)}$) and experimental ($f_{exp}^{(i)}$) spectra are given by:

$$D_1^{(i)} = \frac{1}{2} \int \left| \frac{1}{c_{th}^i} f_{th}^{(i)}(e) - \frac{1}{c_{exp}^i} f_{exp}^{(i)}(e) \right| de,$$

$$D_2^{(i)} = \frac{1}{2} \int \left| \frac{1}{C_{th}^i} F_{th}^{(i)}(e) - \frac{1}{C_{exp}^i} F_{exp}^{(i)}(e) \right| de,$$

with each time the normalization factors: $c = \int_{E_{min}}^{E_{max}} f(e) de$

The F are the integrals of f . The R factor is a conventional one (see E. Zanazzi and F. Jona, Surf. Sci. 62, 61 (1977)) given by:

$$R_X^{(i)} = \frac{13}{6} \frac{e \sum \left| c_i f_{th}^{(i)}(e) - f_{exp}^{(i)}(e) \right|^2}{\sum_e \left| f_{exp}^{(i)}(e) \right|^2}, \quad \text{with } c_i \text{ such that: } \frac{\partial R_X^{(i)}}{\partial c_i} = 0$$

Then for n different spectra, one gets the total metric distance or R factor by:

$$D_{1,2} = \sum_{i=1,n} p_i D_{1,2}^{(i)} \quad \text{and} \quad R_X = \frac{6}{13} \left(\frac{3}{2n} + \frac{2}{3} \right) \sum_{i=1,n} p_i R_X^{(i)}$$

In which p_i is the relative weight for the spectra "i" given by $p_i = \frac{E_{max}^{(i)} - E_{min}^{(i)}}{\sum_{j=1,n} E_{max}^{(j)} - E_{min}^{(j)}}$

It is possible to vary some parameters in order to optimize the agreement between calculation and experiment. Then results are given for a multi-dimensional grid of parameters. The values of the metric distances are given in a special output file. By default for a complete calculation, the file name is the fdmnes conventional output file name with the suffix *_fit.txt*. For a calculation starting calculating the convoluted spectra, the default output file name is *fdmfit_out.txt*. This name can be modified using a specific keyword.

To do this task different keyword must be added in the main indata file:

Experiment <i>Nom_exp</i>	→ Keyword preceding the file names containing the experimental spectra.
Gen_shift <i>7108 7114 21</i>	→ Minimum and maximum energy shift between calculation and experiment and number of value to test
Metric_out <i>File_name_fit.txt</i>	→ Keyword preceding the output file name containing the metric distances
Parameter <i>Par_Gamma_max</i>	→ Keyword preceding any group of correlated parameters → Keyword to specify the convolution broadening Γ_m as a parameter

15. 25. 5 → First and last values and number of values

Parameter

Par_Posx

-0.05 0.05 3

Par_Posy

-0.05 0.05 3

In the previous example the parameters *Posx* and *Posy* of the atom 1 are completely correlated. Thus the atom is displaced along the diagonal.

The parameters can be fitted are:

For the convolution:

<i>Par_ecent</i>	→ Central energy for the arctangent
<i>Par_elarg</i>	→ Energy width for the arctangent
<i>Par_efermi</i>	→ Fermi (or cutting) energy
<i>Par_gamma_hole</i>	→ Hole width
<i>Par_gamma_max</i>	→ Maximum width for the final states
<i>Par_gauss</i>	→ Gaussian width (or resolution)
<i>Par_shift</i>	→ Energy shift
<i>Par_aseah</i>	→ First parameter of the Seah-Dench formula

For the spectra calculation:

<i>Par_a</i>	→ Contraction or expansion of the mesh parameter a in %
<i>Par_b</i>	→ Contraction or expansion of the mesh parameter b in %
<i>Par_c</i>	→ Contraction or expansion of the mesh parameter c in %
<i>Par_abc</i>	→ General contraction or expansion in %
<i>Par_anga</i>	→ Value of the unit mesh angle α
<i>Par_angb</i>	→ Value of the unit mesh angle β
<i>Par_angc</i>	→ Value of the unit mesh angle γ
<i>Par_poporb</i>	→ Orbital occupancy
<i>Par_posx</i>	→ Atom position along x
<i>Par_posy</i>	→ Atom position along y
<i>Par_posz</i>	→ Atom position along z
<i>Par_occup</i>	→ occupancy of the atom
<i>Par_dposx</i>	→ shift of the atom position along x from the original position
<i>Par_dposy</i>	→ shift of the atom position along y from the original position
<i>Par_dposz</i>	→ shift of the atom position along z from the original position
<i>Par_theta</i>	→ position along θ for an atom in spherical coordinate
<i>Par_phi</i>	→ position along ϕ for an atom in spherical or cylindrical coordinate

Under each parameter must be written the first and last values of the parameter followed by the number of value. For the parameters *Par_posx*, *Par_posy*, *Par_posz*, *Par_theta*, *Par_phi* The number of the atom must also be specified. Under the parameter *Par_poporb*, *Par_occup* must be also set the number of the orbital found in the order given under the keyword *atom*. For example, to specify an exchange of charge between two atoms (titanium and oxygen), one has to write:

Atom

22 3 2 2. → 3d titanium orbital with 2 electrons for the initial occupancy
8 2 1 4. → 2p oxygen orbital with 4 electrons for the initial occupancy

Parameter**Par_poporb**

2. 0. 3 1 → The first orbital under *atom* is the 3d titanium.

Par_poporb

4. 6. 3 2 → The second orbital under *atom* is the 4p oxygen.

By default the metric distances are calculated in all the energy range is the intersection between the experimental and calculated spectra. It is possible to cut the lower or and the higher energy part of the spectra by the use of the keyword:

Emin

-10. → Minimum energy for all the spectra

Emax

100. → Maximum energy for all the spectra

It is possible to have different values for the different spectra:

Emin

-10. -5. -20. -20. → Minimum energy for each spectra

Emax

45. 100. 100. 100. → Maximum energy for each spectra

If the energy of the experimental spectra is in keV and not in eV, put the keyword:

Kev

When there are several spectra to compare simultaneously, for the calculation they must be in the same output file. One then must give the number of the column contains the calculated spectra and associate it at the corresponding experimental file. This is done adding a line after each experimental file contains the number of the column. If there are 2 numbers, the first one is the number of the column in the experimental file, the second one being the number of the column in the calculated file. When there is no number, this means than in both experimental and calculated file, the spectra are in the second column, the first one being the energy. Example:

Experiment → keyword

Nom_exp_1.txt → Name of the file containing the first experimental spectra

2 → Number of the column in the calculated file containing the corresponding spectra

Nom_exp_2.txt → Name of the file containing the second experimental spectra

3 → Number of the column in the calculated file containing the corresponding spectra

Nom_exp_2.txt → Name of the file containing the third experimental spectra

4 3 → Numbers of the column in the experimental file *Nom_exp_2.txt*, followed by the number of the column in the calculated file

By default, three confidence factors are used. A fourth one “Rxg” can also be calculated. It is equivalent to Rx but with a unique normalization factor between experiment and calculation. When there is a single spectra, Rxg is equivalent to Rx. To get this confidence factor put the keyword:

Rxg

It is also possible to compare without « fit » the spectra, which can be coming from 2 calculations. In that case, they must have the same number of columns, the first one being the energy. The spectra from the same columns in both files are compared. For this purpose, one must use the keyword « *File_met* » followed by the 2 file names:

File_met

File_1_conv → name of file number 1

File_2_conv → name of file number 2

The keyword *Gen_shift* or *Metric_out* can be used with *File_met* as usual. This keyword can be used alone when both files are existing.

E- Extraction of DAFS scans and spectra

When realizing a DAFS simulation with azimuth dependence (scan), the output scan file after convolution contains the intensity of all the reflections at all energies and all azimuth angles. Often, it is useful to extract from this big file some spectra at specific azimuth angle or some scan at specific energy. It is what is done in this part with the following indata file:

!Indata file for FDMNES, Selection part

Selec_inp → keyword for the input file name (output of the convolution part)
xanout/fe3o4/2008_bland/cc_1221_1221_orig_1m10_scan_conv.txt

Selec_out → keyword for the output file name
xanout/test/fe3o4_scan_selec_conv.txt

Energy → keyword for the selected energy for the scan
 4. → value of the energy for the scan in eV.

Reflection → keyword for the selected reflections
 2 5 6 9 → number of the selected reflections

End

The azimuth scan is extracted at the energy (in eV) given after keyword “*Energy*”.

It is also possible to extract spectra at a given azimuth angle. For this, no keyword “*Energy*” but the keyword “*Azimuth*” and next line the value of the selected azimuth:

Azimuth
 30. → Value of the azimuth

It is possible to have several reflections in the output file at one energy (or one azimuth). It is also possible to have several energies but only one reflection for the scan or several azimuth angles but one reflection for the spectra.

Energy
 4 7.5 → selected energy for the scan

Reflection
 2 → reflection number 2 selected

or:

Reflection
 2 → reflection number 2 selected

Azimuth
 30. 60. 90. → Azimuth selected

Without the keyword « *Reflection* », all reflections are given when there is only one energy or one azimuth asked.

F – Unit cell multiplication

FDMNES contains a tool to build the list of atoms of a unit cell superstructure from the atom list of a simple unit cell. This can be used as the starting list corresponding to a unit cell corresponding to a decrease in symmetry or to specific magnetic order.

Example:

Filout → Output file in which is given the list of atom of the superstructure unit cell
xanout/test_stand/mult_out

Mult_cell → Multiplication factor along a, b and c
 2 2 1 → Here one builds a (2a, 2b, c) unit cell

Atomic_nu → Atomic number corresponding to the atom type under « *Unit_cell* »
 26 38 33

Unit_cell → to introduce the original unit cell
 3.92430 3.92430 12.36440 90. 90. 90. → a, b, c, α , β , γ
 1 0.50000 0.00000 0.25000 → list of all the atoms of the unit cell
 1 0.00000 0.50000 0.25000
 1 0.50000 0.00000 0.75000
 1 0.00000 0.50000 0.75000
 2 0.00000 0.00000 0.00000
 2 0.50000 0.50000 0.50000
 3 0.00000 0.00000 0.36000
 3 0.00000 0.00000 0.64000
 3 0.50000 0.50000 0.86000
 3 0.50000 0.50000 0.14000

End

The original unit cell is introduced with the keyword « *Unit_cell* » under which are given the unit cell parameters, and the list and the list of all the atoms in it. In first column is set the atomic number, but if using the optional keyword « *Atomic_nu* ». In this case, first column contains the atom type number and the correspondence between atom type number and atomic number is given under « *Atomic_nu* ».

G - FDMX User's guide

FDMX is an extension from J. Bourke and C. Chantler, University of Melbourne, Australia. When using it, thanks to cite:

Jay Daniel Bourke, Christopher Thomas Chantler and Yves Joly
"Extended X-ray Absorption Fine Structure Calculations Using the Finite Difference Method"
 J. Synchrotron Rad. **23**, 551-559 (2016).

FDMX is an enhanced approach to calculating both XANES and XAFS spectra using the finite difference method and the core routines of FDMNES. The easiest way to use FDMX is to simply include the keyword ***fdmx*** in your regular input file, and the code will automatically optimize computational parameters to ensure accuracy over a wide energy range. FDMX will process the spectrum to include thermal effects and electron/hole lifetimes, and add background absorption from more loosely bound electrons. The code is currently designed for use with K-edge calculations only, however most functionality will also work with other edges (note keywords for hole widths and background absorption). Additional (optional) controls for FDMX are possible with specific keywords.

As calculations using FDMX may take several hours, particularly for structures with few or no axes of symmetry, it is strongly recommended that you compile and run the code using the MUMPS libraries.

fdmx

Activates FDMX computation and optimization of parameters with respect to energy, allowing for accurate EXAFS spectra, and triggers the processing routine at the end of the calculation. When this keyword is used, *the Radius keyword is no longer required*.

fdmx_proc

Use to activate a post-processing routine from FDMX, implementing thermal, inelastic scattering, background, and other effects without explicitly calculating the absorption spectra. This requires an existing output file with absorption cross sections from a previous calculation.

E_cut

4.0 ! Val

Similar to the keyword for ***convolution*** in FDMNES, this will suppress transitions below the nominated Fermi level, expressed in eV. Some absorption will still be possible below this level due to the hole width and IMFP effects, and background absorption from other edges. **Use of this keyword is strongly recommended.**

Radius

8.3 100. 7. 250. 5. 400. 4. ! Val, En, Val ...

Allows the user to specify how the cluster radius changes with energy. By default when using *fdmx* procedure, the radius is set (as above) to 8.3 Å for energies up to 100 eV, then 7 Å for

energies up to 250 eV, then 5 Å etcetera. Using a high or constant radius at high energies may result in very long calculations, while low radii will cause XANES and XAFS structure to be missing.

adimp

0.24 0. 100. 0.20 250. 0.16 400. 0.12 500. 0.08 ! *Val, En, Val ...*

Allows the user to specify how the inter-point distance (grid density) changes with energy. By default, the inter-point distance is set (as above) to 0.25 Å for energies up to 0 eV, then 0.22 Å for energies up to 10 eV, then 0.18 Å etcetera. Using a high or constant inter-point distance at high energies may produce inaccurate results, while low values will lead to long calculations.

Gamma_hole

2.3 ! *Val*

By default, FDMX will include a core-hole relaxation based on K-shell tabulations from Scofield and Kostroun (Z=21-50) and “B” (Z=51-100)? Alternatively one may provide an explicit core-hole lifetime (in eV) with the keyword ***Gamma_hole*** or suppress this effect by using the keyword:

nohole

It is required that you use either ***Gamma_hole*** or ***nohole*** for calculations involving edges other than K.

IMFPin

imfpdatafile.txt ! *Filename*

By default, FDMX will include photoelectron lifetime broadening based on the electron inelastic mean free path (IMFP) of the absorbing material. For elements the IMFP values are tabulated from optical data model calculations, while for compounds they are estimated using the TPP-2M equation. Users may instead include their own IMFP values with the ***IMFPin*** keyword followed by a filename that contains IMFPs in the format:

En Val

where *En* is in eV and *Val* is in Å. Alternatively, IMFPs can be calculated directly from optical energy loss data (i.e. optical ELF), which may be obtained from band structure packages such as WIEN2k (in WIEN2k this data is obtainable using the *optic* routines and is output to a file named *case.eloss*). To use ELF data to compute the IMFP, use the keyword:

ELFin

elfdatafile.txt

where the data is formatted the same as for an IMFP input file, excepting that the values are unitless. The use of ***ELFin*** will overwrite other IMFP values up to the maximum energy of the ELF provided (or 120 eV, whichever is lower), after which default IMFP values will be used. By default, the IMFP is calculated from the ELF using a Lindhard (fast) representation of the

dielectric function. However, it is recommended that you use an additional keyword for a Mermin representation:

Mermin

1 ! Val

The number following the ***Mermin*** keyword specifies an iteration value, with 0 corresponding to a Lindhard representation, 1 to a standard Mermin representation (recommended), and higher values to successive iterations of a self-consistent Mermin model. This calculation may take some hours, and iteration values higher than 3 will not normally be useful for XAFS. For more information see Bourke & Chantler, *J. Phys. Chem. Lett.* **6** 314 (2015).

Use of the ***ELFin*** keyword will generate an additional output file with the extension _ELF.txt containing the energy- and momentum-dependent ELF of the system. Use of ***Mermin*** without ***ELFin*** will produce no effect.

noIMFP

Use to omit the effects of photoelectron lifetime, and hence represent an infinite IMFP.

Tmeas

298 ! Val

Use to define the temperature, in Kelvin, at which the absorption is being calculated. This value will affect the Debye-Waller Factor, and is set to 0 by default.

TDebye

345 ! Val

Use to define a Debye temperature (in Kelvin). This is strongly recommended for compound structures.

DWfactor

0.13 ! Val

Use to define an effective Debye-Waller factor (in Å). By default, this will otherwise be calculated based on the dominant oscillators contributing to the high energy XAFS spectrum.

noDW

Use to omit thermal effects – effectively a zero Debye-Waller factor.

noBG

By default, FDMX will include background absorption from more loosely bound electrons (i.e. L, M edges etc). Use this keyword to omit such contributions. Note that the background absorption currently only works properly for K-edge calculations, and therefore for other edges this keyword is required.

For comparison with experimental data, it is useful to define an extra function to address the edge jump discrepancy (triangle effect). This function comes in two flavours, and is implemented by the following relation:

$$\mu_0(E) \rightarrow \mu_0(E)(1 + f(E))$$

where μ_0 is the atomic (or atom-like) background absorption.

Victoreen

0.8 0.4 ! Val A, Val B

Adds a background function from the Victoreen Equation with A and B parameters following:

$$f(E) = A \left(\frac{E_0^3}{E^3} \right) + B \left(\frac{E_0^4}{E^4} \right)$$

This applies to all energies above the absorption edge E_0 . $A+B$ provides the relative edge jump correction, while A/B controls the rate of decay.

Expntl

0.1 100. ! Val A, Val B

Adds a background function of an exponential form with A and B parameters following:

$$f(E) = Ae^{-\frac{(E-E_0)}{B}}$$

This applies to all energies above the absorption edge E_0 . A provides the relative edge jump correction, while B controls the rate of decay.

cm2g

Output results in units of cm^2/g instead of Mbarn.

H - 2D Diffraction

FDMNES can simulate 2D – resonant diffraction experiments. The corresponding output are (h,k,l) peaks with non-integer indices. Truncation rods at a specific energy can be given as well as spectra at specific points of the truncation rods.

1) *Surface structure description*

The surface is supposed perpendicular to the (a,b) plane, thus eventually, the unit cell must be especially written in this way. Note that c is not necessarily perpendicular to the surface. The surface to simulate must contain some atomic layers that we call “surface” and most often a substrate taken as semi-infinite crystal that we call “bulk”.

The substrate is taken into account with the keyword:

```
Bulk                                ! Structure of the substrate
      8.1724 8.1724 4.0862 90. 90. 90. ! a, b, c, alpha, beta, gamma
47 0.0 0.0 0.0                       ! Z, x, y, z
47 0.5 0.5 0.0
47 0.5 0.0 0.5
47 0.0 0.5 0.5
```

On top of the bulk atoms, the surface atoms are set using keyword “*Surface*”:

```
Surface
      8.1724 8.1724 4.0862 90. 90. 90. = a, b, c, alpha, beta, gamma
27 0.0 0.5 0.2                       ! Z, x, y, z
.....
```

Note that the unit cell parameters are not necessary the same than the bulk unit cell parameters. The distances between the top most bulk atom and the bottom less atom of the surface is given by the lowest z value (in the 4th column). This can be change using keyword “*Surface_shift*”. When using it, this distance is given under this keyword:

```
Surface_shift
1.25          → distance between bottom most surface atom and top most bulk layer (A)
```

To give a shift in the 3 directions one writes:

```
Surface_shift
0 0.23 1.25    → shifts along a and b and distance between bulk and surface (A)
```

When there is many surface atoms, it can be convenient to consider it as a film which is considered as “bulk-like” along some unit cell length. In such a film the atomic form factor are calculated as in a 3D material. For this purpose one use the keyword « *Film* »:

Film → Set a simulation on a 2D film, followed by the unit cell description
 8.1724 8.1724 4.0862 90. 90. 90. = *a, b, c, alpha, beta, gamma*
 47 0.0 0.0 0.0 ! Z, x, y, z

By default the film thickness is equal to the '*c*' unit cell parameter, eventually multiplied by the cosines of the angle between *c* and the basal plane. One can use a specific thickness using the « *Film_thickness* » keyword, followed by its value in Angstrom not necessarily equal to a multiple of *c*.

Film_thickness
 75

The interlayer distance between the topmost substrate layer and the bottom film layer is by default zero. This can (and must) be changed with the following keyword:

Film_shift
 3.5 → *z* when only one number is given (in Angström), *z* is perpendicular to the surface, direction outwards.

Film_shift
 0. 0. 3.5 → *x, y, z*, all in Angström, *x* and *y* are along *a* and *b* of the film.

Film_shift
 0. 0. 3.5 45. ! *x, y, z, angle*

A fourth number can be added. It is the angle between the *a* film axis and the *a* bulk axis. Can be useful for epitaxial growth when the a film corresponds to the (*a,b*) diagonal of the bulk

By default, the bottom atom in the film is the one defined with the lowest *z* position under the keyword "Crystal". This can be changed defining a new zero, in reduced unit. The atoms between *z* = 0 and this value are just shifted by +1.

Film_zero
 0.25 → The bottom atom will be the next one with *z* >= 0.25. Then the stacking is as usual.

Note that the keyword "*Spgroup*" can be used for a film, but not for the description of "*Surface*" or "*Interface*" as seen bellow. It is possible to use only a film and a bulk. To improve the calculation it is possible to add some surface layer atoms on top of the film. For this purpose, we just use the same "*Surface*" keyword as before. Because there is a "film", the code understands that the "surface" is on top of the film and not on top of the bulk.

It is also possible to consider an interface between the bulk and the film. This is done with the same format than with the surface, film or bulk atoms but with the keyword:

Interface
 8.1724 8.1724 4.122 90. 90. 90. = *a, b, c, alpha, beta, gamma*
 26 0.0 0.0 0.0

The interface distance from the bulk is given in the same format than with the “surface” using the keyword:

Inter_shift

1.26

One can add a cap layer with:

Cap_layer

2.8837228875 2.8837228875 7.0636496 90. 90. 120.
79 0. 0. 0. 1.
79 0.6666666667 0.3333333333 0.3333333333
79 0.3333333333 0.6666666667 0.6666666667

Cap_thickness

17.1

A Debye disorder in Angström can be given using « *Cap_B_iso* » or « *Cap_U_iso* ». This disorder is the same for all the atoms of the cap layer:

Cap_B_iso

80.0 $\rightarrow = 8\pi^2 \langle u^2 \rangle$ in Å²

Cap_U_iso

1.0 $\rightarrow = \langle u^2 \rangle$ in Å²

By default the distance between film and cap is the sum of the atom radii. This can be changed with:

Cap_shift

3.0

The roughness of the « bulk » or the interface bulk/film, of the film or the interface film/cap layer and of the cap layer can be taken into account. They are described by the function $0.5 \operatorname{erfc}(z - z_0/\sqrt{2}R)$, where z_0 is the origin taken at the middle in the inter-layer above the last atomic layer. This is done with :

Bulk_roughness

1.5

\rightarrow roughness at the bulk surface or at the interface bulk/film

$\rightarrow R$

Film_roughness

4.

\rightarrow roughness at the film surface or at the interface film/cap layer

Cap_roughness

17.1

\rightarrow roughness at the cap layer surface

With the keywords "*Occupancy*" and/or "*Atom_B_iso*" (or "*Atom_U_iso*") another 1 or 2 other columns, in the 5th and 6th columns, can be added for surface, film, inter-face and bulk slabs, containing the occupancy rate (between 0 and 1) and the Debye parameter, $8\pi^2\langle u^2 \rangle$, in Å². The column order of these parameter corresponds to the order of these keywords in the indata file.


Atom_B_iso

```
47 0.0 0.0 0.0 0.82      ! Z, x, y, z
....
```

Occupancy

Atom_B_iso

```
47 0.0 0.0 0.0 0.8 0.82      ! Z, x, y, z, Occupancy, Debye
```

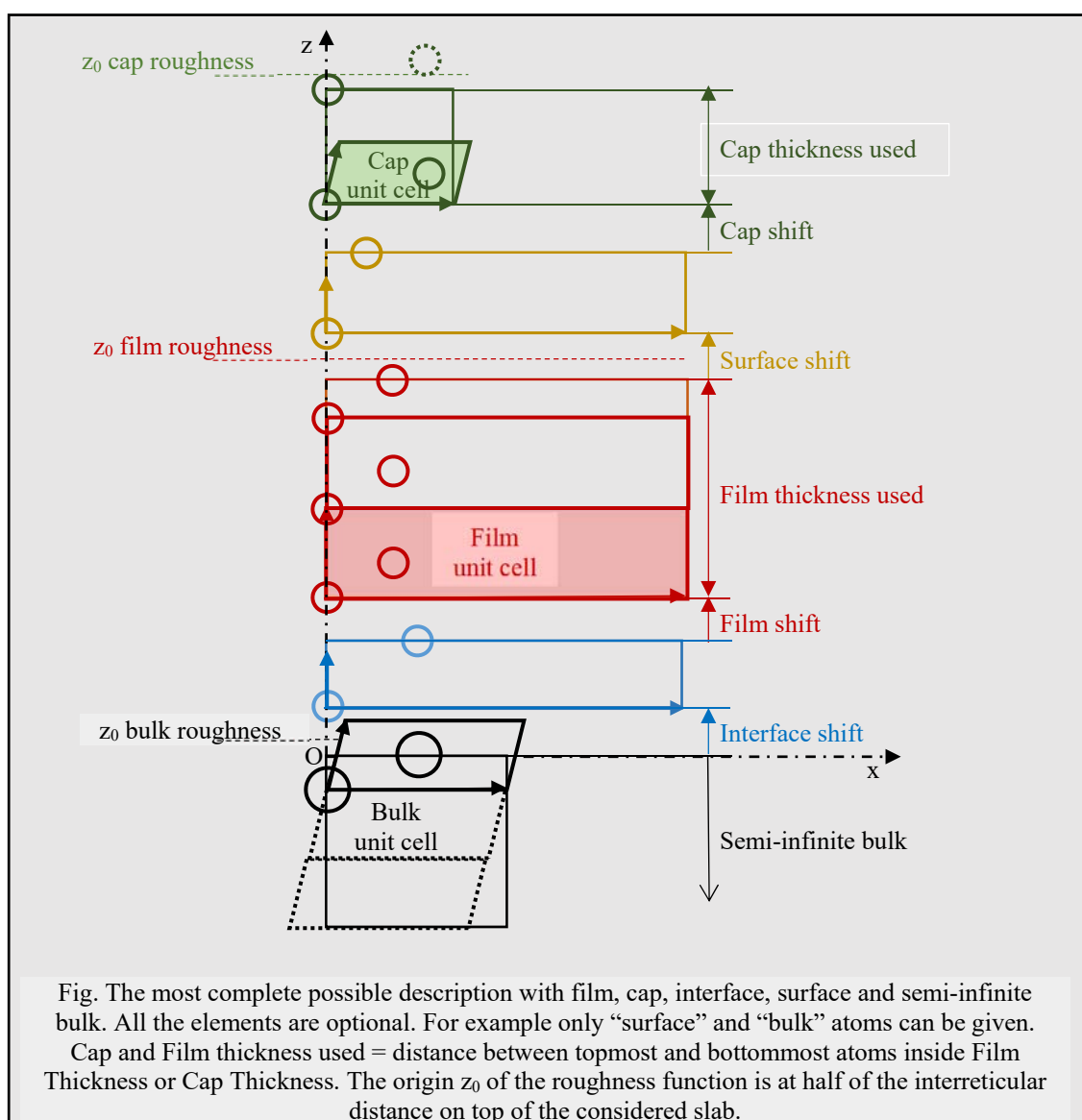

2) Reflection choice

DAFS

```
0 0 0.04 0.04 8.16 1 1 0.      ! h, k, l_min, l_step, l_max, sigma, sigma, azimuth
0 0 8.18 0.002 8.24 1 1 0.
0 0 8.28 0.04 10.00 1 1 0.
```

There is thus a scan along the rod from l_{\min} of to l_{\max} .

Note that h and k can be used as non-integer when the (a,b) film periodicity is rational (and not equal) to the bulk periodicity.

It can be more convenient to use the experimental operation mode as the one defined in « Spec ». In that case use :

DAFS_2D

```
0 0 0. 0.05 5. 4 1 4 0 0 0. 0.      ! h, k, l_min, l_step, l_max, Operation mode, Angle
0 0 0. 0.05 5. 4 1 4 0 0 90. 0.      ! h, k, l_min, l_step, l_max, Operation mode, Angle
0 1 0. 0.05 5. 2 1 2 0 0 30. 10.
0 0 1.          4 1 4 0 0 90. 0.      ! h, k, l, Operation mode, Angle
```

After h,k,l (with or without scan along a truncature rod), one gives the 5 integer of the operation mode of the same table than in « Spec » (and nearly corresponding to the You'th paper (J. Appl. Cryst, 32, 614 (1999)), then the values of the fixed angles :

	Detector	Reference	Sample	Sample	Sample
1	delta	alfa = beta	eta	eta	eta
2	nu	alfa	mu	mu	mu
3	Qaz	beta	chi	chi	chi
4	Naz	Psi	phi	phi	phi
5	x	X	eta = delta / 2	eta = delta / 2	eta = delta / 2
6	x	X	mu = nu / 2	mu = nu / 2	mu = nu / 2

delta = δ : detector polar angle
 nu = ν : azimuthal detector angle (independent from mu)
 Qaz : azimuth of the rotation of the diffraction vector / laboratory
 Naz : azimuth of the rotation of the normal to the surface / laboratory. Naz = 0 when normal is horizontal. Specular case: polarizations σ et $\pi \rightarrow$ Naz = 90 and 0°.
 alfa = α : incidence angle
 beta = β : outgoing angle
 psi = ψ : azimuthal angle of the sample rotation around the diffraction vector. Its zero is when the reference vector is in the incidence plane in the positive direction of the incoming beam; diffraction vector is then horizontal. By default the reference vector (in reciprocal space) is (0,0,1). To change it, see keyword "Setaz".
 eta = η : angle polaire de l'échantillon
 mu = μ : angle azimuthal de l'échantillon

chi = χ : autre angle de l'échantillon (chi = 0, correspond au « z-axis mode »)
 phi = ϕ : angle azimutal de l'échantillon

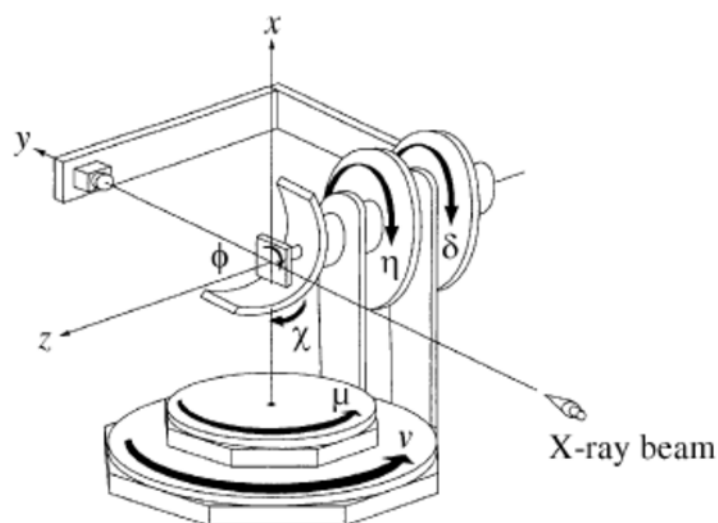


Figure from You's paper showing the angle definitions corresponding to a « 4+2 » diffractometer when all angles are 0

For example, the operation mode *4 1 4 0 0* means « Naz fixed, alfa=beta, phi fixed ».

The following angles are the values of the fixed angles. For example :

« *4 1 4 0 0 0. 0.* » means « Naz =0., alfa=beta, phi =0 ».

« *4 1 4 0 0 90. 0.* » means « Naz =0., alfa=beta, phi =90 ».

« *2 1 2 0 0 30. 10.* » means « nu =30., alfa=beta, mu =90 ».

Some operation modes are not possible or not coded :

<i>4 0 i j 0</i>	for any i and j
<i>i 0 1 3 0</i>	for i = 1, 2 or 3
<i>i 0 1 4 0</i>	for i = 1, 2 or 3

For the specular reflection, one may only use :

<i>4 1 4 0 0</i>	
<i>0 1 3 4 0</i>	with only chi = 0 or 90
<i>0 1 4 3 0</i>	with only chi = 0 or 90

When the experimental orientation matrix (U) times the (h,k,l)→Q matrix transformation (B) is known (the so call UB matrix), this one can be used:

Mat_UB

```
-0.6713131898 -0.331076068 -0.005363453371
0.331087442 -0.6713288502 -0.0004569418448
-0.004608158701 -0.002782139698 0.7485135622
```

By default, the reference axis from which the psi angle is defined is the normal to the surface. It is possible to define another direction, as in the experiment using:

Setaz

1 1 0 → direction (in reciprocal space coordinate)

It can be interesting in the convolution part of the calculation to get not only the spectra but also the scan along the rod at specific energies. To get the corresponding new file (with extension "_tr"), write:

Transpose_file

15. 17.

The new file is in fact just a transpose of the standard output file, taken in the example at the energies 15 and 17 eV.

The film, surface, interface and bulk unit cell parameters can be all different. By default the (h,k,l) indices correspond to the bulk parameters. One can make that it corresponds to the film ones with keyword:

hkl_film

Note that when film and unit cell parameters are not rational between them, film and bulk non specular reflections do not interfere.

By default the code calculate the intensity for well-defined outgoing polarization. Using “*DAFS_2D*” keyword, one thus get output spectra for both σ and π outgoing conditions. To get only the sum of them, that is as being without analyzer, use the keyword:

No_analyzer

In SRXRD, it is possible to consider a Helmholtz layer. See the keyword “*Helmholtz*” in another section of this manual.

Because the bulk is often a dense simple metal, it can be faster to simplify the calculation in this part. It is possible to avoid the SCF calculation and/or the finite difference method. For this use:

No_SCF_bulk

→ no SCF calculation in the bulk, whatever it is SCF or not in the surface

Green_bulk

→ in the bulk, the calculation is done using the multiple scattering theory, whatever is the calculation mode in the surface.

List of *fdmnes* keywords

Abs_B_iso	58	Convolution	53, 55
Abs_U_iso	58	COOP_atom	44
Abs_before	54	COOP_dist	44
Absorber	17, 39, 42	Core_resolved	51
Adimp	48, 49, 72	Crystal	13, 27, 34
All_conv	54	Crystal_p	15
All_nrixs	21	Dafs	22, 23, 25, 79
Allsite	41	DAFS_2D	79
Ang_spin	29	Dafs_exp	24
Atom	15, 27, 28, 34, 64	Dafs_exp_type	24
Atom_B_iso	26, 78	Dead_layer	60
Atom_conf	16, 35	Density	38
Atom_U_iso	26, 78	Density_comp	39
Atom_nsph	35	Density_all	39
Atomic_nu	69	Dec	55
Axe_spin	29	Delta_E_conv	30
Azimuth	67	Dilatorb	36
Basecomp	46	Dip_rel	19
Bulk	75	Dipmag	18
Bulk_roughness	77	Directory	54
Calculation	53	Doping	20
Cap_B_iso	77	Double_cor	60
Cap_layer	77	Dpos	40
Cap_roughness	77	DWfactor	73
Cap_shift	77	E_cut	54
Cap_thickness	77	E_out_max	46
Cap_U_iso	77	E1E2	18
Cartesian	40	E1E3	18
Center	38	E1M1	18
Center_abs	38	E2E2	18
Center_s	38	E3E3	18
Check	45	Ecent	56
Check_all	45	Edge	18
Check_biref	61	Eimag	36
Check_coabs	45	Elarg	56
Check_conv	56	Emax	65
Check_mat	45	Emin	65
Check_pot	45	End	11
Check_sph	45	Endjump	43
Check_tens	45	Energpho	18
Chlibre	35	Energy	67
Circular	61	Epsii	41
Clementi	34	Estart	57
cm2g	74	ELFin	72
Comment	43	Excited	48
Conv_gaus	62	Experiment	63, 65
Conv_out	50	Expntl	74

Extract	20	M1M1	18
Extract_tens	40	Magnetism	26, 42
Extractpos	41	Mat_polar	20
Extractsym	41	Mat_UB	80
File_in	62	Memory_save	51
File_out, Filout	13, 62, 69	Mermin	73
File_met	66	Metric_out	63
Fdmx	71	Molecule	13, 32
Fdmx_proc	71	Muffintin	48
Film	14, 76	Mult_cell	69
Film_Cif_file	15	N_self	30
Film_Pdb_file	15	No_analyzer	81
Film_roughness	77	No_check	44
Film_shift	76	No_E1E1	18
Film_thickness	76	No_E1E2	18
Film_zero	76	No_E1E3	18
Flapw	42	No_E2E2	18
Flapw_n	42	No_Fermi	29
Flapw_n_p	43	No_res_mag	25
Flapw_psi	41	No_res_mom	25
Forbidden	59	No_SCF_bulk	81
Fprime	60	noBG	74
Fprime_atom	60	noDW	73
Full_self_abs	25	noIMFP	73
Gamma_hole	56	Noncentre	38
Gamma_var	56	Nonexc	48
Gamma_max	56	Nonrelat	27
Gaussian	56	Normaltau	48
Gen_shift	65	Norman	50
Green	18	Nrato	43
Green_bulk	81	NRIXS	21
Hedin	34, 41	NRIXS_mono	21
Helmholtz	32	Occupancy	15, 26, 78
hkl_film	81	Octupole	18
Hubbard	31	One_run	44
IMFPin	72	Overad	40
Interface	76	Overlap	50
Inter_shift	77	P_self	30
iord	48	P_self_max	30
Jump	43	Parameter	63
Kev	65	Perdew	32
Length_line	51	Photoemission	55
Lmax	49	Polarize	19
Lmax_nrixs	21	Quadrupole	18
Lmax_tddft	31	R_self	30
Lmaxfree	49	Range	17
Lmaxso	10	Rangel	17
Lmaxso_max	49	Raydem	50
Lminus1	37	Radius	13
Lplus1	37	Rchimp	51

Reflection	67	Spherical	39
Relativism	27	Spinorbit	26
Rmtg_Z	50	Step_azim	25
Rmt	48	Stokes	61
Rmtv0	50	Stokes_name	61
Rotsup	41	Supermuf	50
RPA	31	Surface	75
Rpotmax	37	Surface_plane	60
Rxg	66	Surface_shift	75
Rydberg	37	Sym	45
S0_2	57	Symsite	46
Sample_thickness	58	Table	56
Scan	59	TDDFT	31
Scan_conv	59	TDebye	73
Scan_file	59	Thomson	59
SCF	29	Tmeas	73
SCF_exc	30	Trace	38
SCF_mag_free	31	Transpose_file	81
Screening	32	Unit_cell	69
Seah	56	V0imp	36
Selec_core	58	Victoreen	74
Selec_inp	67	Vmax	36
Selec_out	67	Xalpha	32
Self_abs	24	Xan_atom	43
Setaz	80	Z_absorber	17
Spgroup	14	Zero_azim	23
Sphere_all	39		